




# Precision performance of Crystal Digital PCR™ vs. quantitative PCR

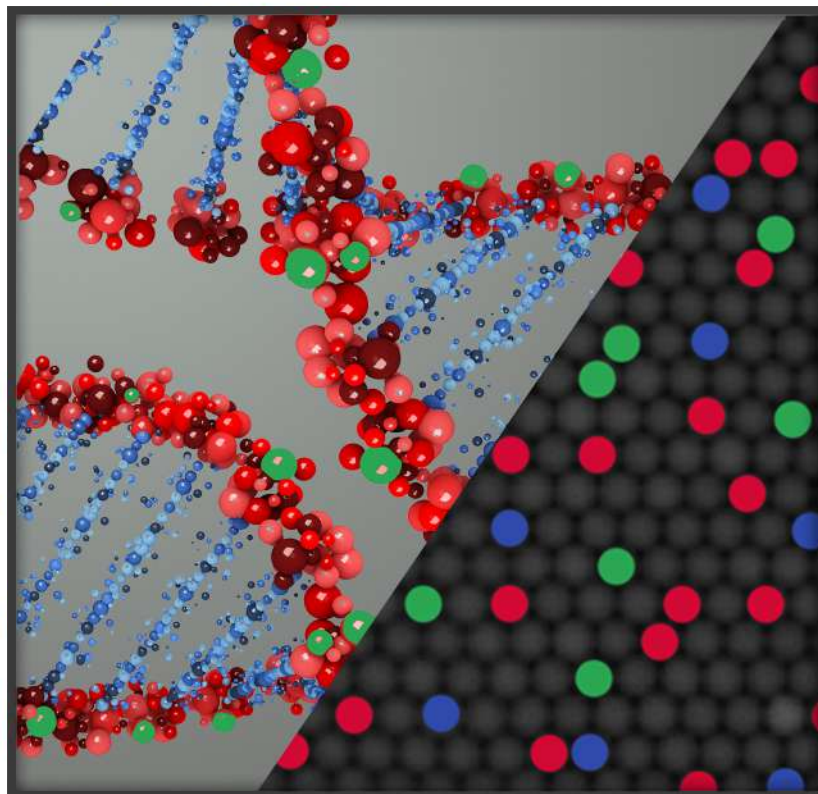

## Experimental Setup

Real-time PCR quantification, also called quantitative PCR (qPCR) compares cycle quantification (Cq) values of unknown samples to that of standards containing known amounts of nucleic acids. Contrary to qPCR, digital PCR (dPCR) enables absolute quantification without the need for a reference standard. Digital PCR systems claim to provide higher precision and reproducibility than qPCR. Here, we compare the precision and reproducibility of the Crystal Digital PCR™ (cdPCR) technology using the Naica™ System to that of qPCR technology. In this study, we made conscious efforts to reduce the experimental variability in order to evaluate only the variability inherent to the technology. We compared the quantification variability of 23 qPCR and cdPCR technical replicates from a single unique PCR master mix spiked with a final concentration of 175 cp/μl of human genomic DNA. The same experimenter and equal final volumes of 25 μL were used to quantify the target gene Albumin (*ALB*) for both techniques (see experimental setup, Fig. 1).

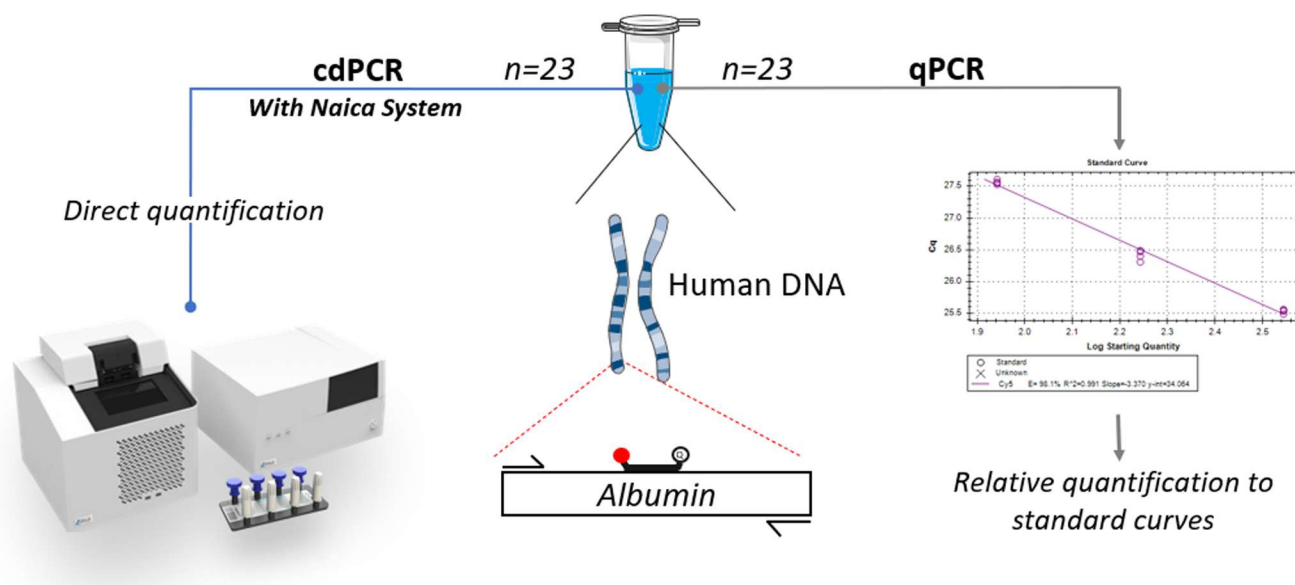


Figure 1: Experiment design to compare the variability of qPCR and Crystal Digital PCR using a customized TaqMan® PCR assay targeting the human Albumin (*ALB*) gene from genomic DNA. A total of 46 data points for *ALB* concentrated at 175 cp/μl was generated from a single PCR master mix, which was divided into 23 replicates quantified in qPCR with standard curves and 23 replicates directly quantified with the Naica System.

## Assay Validation and Quantification

Generally, qPCR standard curves are generated from a different matrix (commercial DNA or standardized samples) than the experimental samples being assayed. The amplification efficiency is not systematically determined for each sample but based instead on a qualitative internal positive control. Thus, the presence of inhibition or non-specific amplification cannot be fully evaluated. In this study, we simultaneously performed standard curves and sample quantification using the same DNA matrix. Standard curve analysis must have a coefficient of determination ( $R^2$ )  $\geq 0.990$  (1). A high efficiency of 98.1% ( $R^2 = 0.991$ ) was obtained (Fig. 2A) within the range of quantification (Fig. 2B).


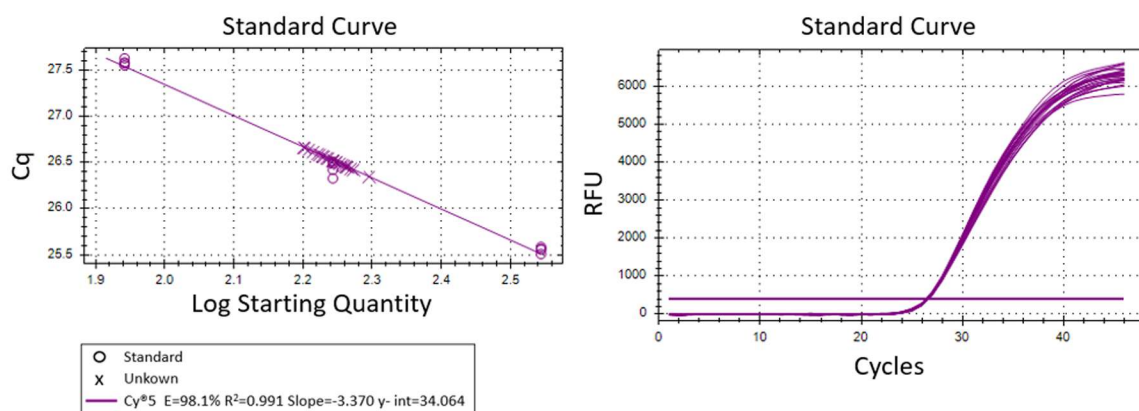


Figure 2: qPCR analysis for 23 data points tested. (A) Standard curve analysis including four replicates per concentration (Efficiency; E).

As cdPCR allows direct quantification, there is no need to utilize standard curves. In addition, cdPCR is less sensitive to PCR inhibiting molecules naturally present in DNA extractions (for more information, see application note about effect of inhibitors on digital PCR at: <https://www.stillatechnologies.com/effect-of-inhibitors-crystal-dpcr-vs-qpcr/>). Indeed, in cdPCR, direct quantification can be determined as long as clear thresholds distinguishing positive and negative clusters can be defined automatically by the Crystal Miner software or manually by the experimenter.


## Comparison of Precision Performance

From the same PCR master mix, 46 technical replicates were created and split into two test sets of 23 sample replicates. These were then measured and compared using qPCR and cdPCR (Fig. 3A). A total of 604,867 droplets were generated for the 23 samples analyzed in cdPCR, thus an average of 26,299 droplets generated per sample. Bessel-corrected coefficient of variation (%CV) was calculated for each set of values (Fig. 3B). The respective CV values show that the measurement variability of cdPCR (%CV= 2.3) is more than 2-fold less than that of qPCR (%CV= 5.0).

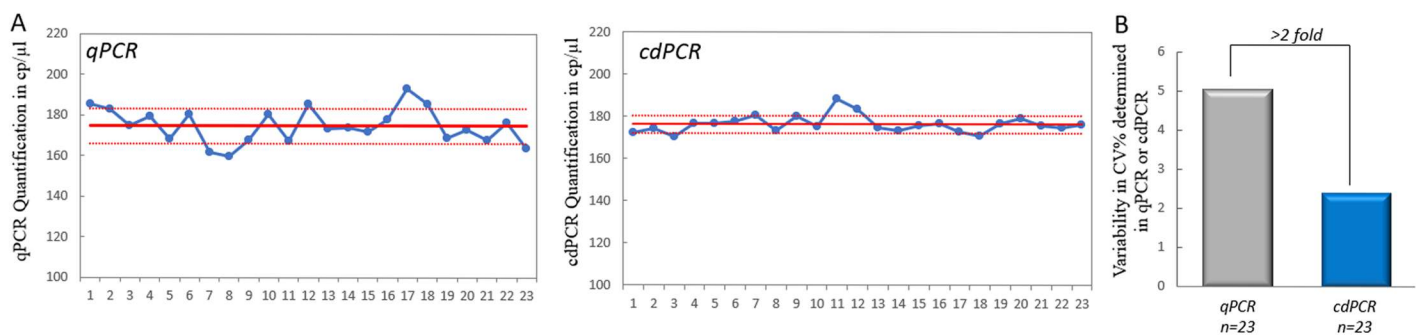




**Figure 3: qPCR and cdPCR comparative analysis. (A)** Graphs showing the quantification of 23 data points generated by qPCR (left) and by cdPCR (right) with one standard deviation (dashed red lines) and the mean quantification value (solid red line) of the 23 replicates. **(B)** Bessel-Corrected CV of the 23 concentration values measured by each technique. CV (%CV=2.3) of cdPCR is 53% less than that of qPCR (%CV=5.0).

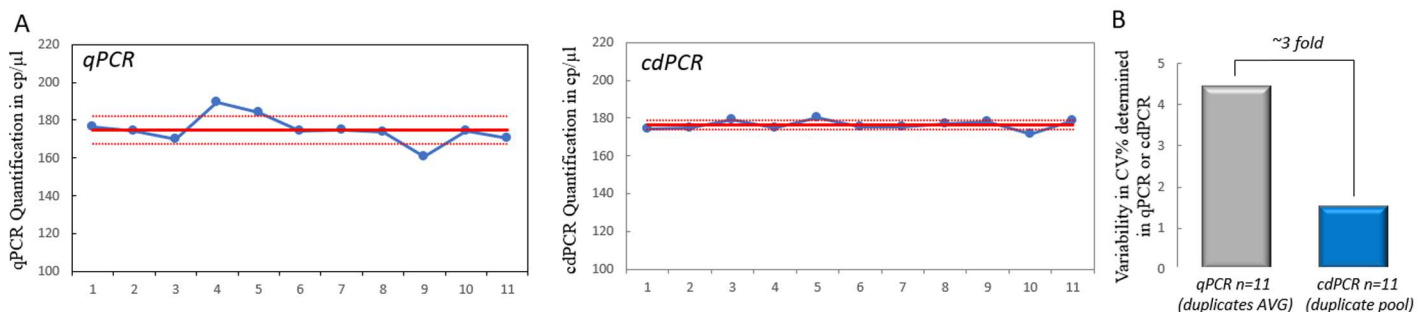Importantly, cdPCR allows sample replicates to be pooled and analyzed as a single larger sample. Here, we compared cdPCR pool in groups of two to the 11 averages of duplicates in qPCR. When cdPCR replicates are pooled, the measurement variability of cdPCR (%CV=1.5) is almost 3-fold less (65.9%) than that of qPCR duplicates average (%CV=4.4) (Fig. 4A and 4B). Indeed, the higher the number of pooled wells, the higher the total analyzed volume and the total number of analyzed droplets, thus lowering both the sampling error and the partitioning error, resulting in a decrease in the quantification uncertainty.



**Figure 4: qPCR and cdPCR comparative analyses where cdPCR replicates were randomly pooled in groups of two and compared to qPCR replicates randomly averaged in groups of two (more information on pooling method at: <https://www.gene-pi.com/item/pooling>).** (A) Graphs showing the average of 11 groups of two (duplicates) quantified using qPCR (left) and the first 11 two-by-two pooled replicates quantified using cdPCR (right) with one standard deviation (red dashes), and the mean quantification value (solid red line). **(B)** Bessel-corrected CV of the concentration values measured using each technique. CV of cdPCR (%CV = 1.5) is 65.9% less than that of qPCR (%CV = 4.4).

## Technical Note Highlights

-  Crystal Digital PCR displays a 2-fold lower measurement variability than qPCR, owing to its end-point determination, its direct quantification, and the high number of crystal droplets generated per sample.
-  Variability in the shown case study can be further decreased to almost 3-fold for Crystal Digital PCR when two wells are grouped using the pooling tool of the Crystal Miner software.
-  Crystal Digital PCR provides the high quantification precision required for various applications including copy number variation, low-fold change and rare mutation detection.

 To learn more about digital PCR, please visit Stilla Technologies  
  Learning Center at [www.gene-pi.com](http://www.gene-pi.com)

### Contact

### STILLA Technologies

Villejuif Biopark  
1 mail du professeur Georges Mathé  
94 800 Villejuif, France

[info@stilla.fr](mailto:info@stilla.fr)  
[www.stillatechnologies.com](http://www.stillatechnologies.com)

+33 9 82 29 50 50



MKT- 00015 Rev A