Limulus Color KY Single Test Wako**Package Insert****【Summary and History】**

Lipopolysaccharide (LPS) also known as endotoxin, is a component of the outer membrane of Gram-negative bacteria. Since endotoxin causes serious symptoms such as fever and shock in the human body, the testing of bacterial endotoxins in pharmaceuticals (especially parenteral drugs) is important.

In 1956, F. B. Bang reported clotting of the horseshoe crab blood due to Gram negative bacteria¹⁾. In 1964, J. Levin and F. B. Bang discovered that clotting of Limulus Amebocyte Lysate (LAL) is caused by endotoxin²⁾. Subsequently, bacterial endotoxin detection using the lysate is widely used as a sensitive and convenient method. This method, known as the Bacterial Endotoxin Test (BET) was listed in the United States Pharmacopoeia in 1980 and in the Japanese Pharmacopoeia in 1988^{3), 4)}. Currently, there are gel-clot, turbidimetric and chromogenic techniques as methods of the bacterial endotoxin test. A bacterial endotoxin test is to be performed by one of these techniques⁴⁾.

Kakinuma *et al.* reported that the lysate reacts not only to endotoxin but also to (1→3)- β -D-glucan (β -glucan)⁵⁾. Iwanaga *et al.* reported that the cascade system in the lysate activated by β -glucan was different from the one activated by endotoxin^{6), 7)}. It was also described that an LAL-reactive material coeluted from cellulose membranes⁸⁾.

A false positive such as from β -glucan or any substance other than endotoxin can present an issue when testing.

This kit contains Limulus Color KY Single Test vials and a Control Standard Endotoxin (CSE) which can specifically, and rapidly detect and quantify endotoxin with high sensitivity. The Limulus Color KY Single Test is a lyophilized lysate including the synthetic chromogenic substrate and buffer constituents. This lysate is a single test reagent that can be measured simply by adding 0.2 mL of sample to the prefilled lyophilized product vial. The Limulus Color KY Single Test is suitable for kinetic analysis by using the Toxinometer[®] Measurement System (Toxinometer[®]).

The CSE is lyophilized LPS which has been extracted and purified from *E. coli* UKT-B by the phenol method^{9), 10)}, including mannitol and glycine as additives. The CSE contains 500 ng of endotoxin and is evaluated in comparison to the Japanese Pharmacopoeia Reference Standard Endotoxin. The potency is determined and printed on the label as a reference.

【Features】

1. The lysate and chromogenic synthetic substrate are formulated together to form the Limulus Color KY Single Test reagent.
2. The reagent can specifically detect endotoxin with high sensitivity without being affected by β -glucan in the sample.
3. The reagent reacts with endotoxin in a wide range of concentrations, some samples may not need to be diluted.
4. Due to its high sensitivity, sample dilution can reduce the chance of interfering factors in the sample.
5. The Toxinometer[®] can quantify endotoxins with the kinetic-chromogenic technique.
6. Limulus Color KY Single Test vials result in no reagent wasted after testing.

[Principle]

The mechanism of this chromogenic reagent in the presence of endotoxin is summarized in Figure 1. The serine protease precursors contained in the lysate are sequentially activated, the clotting enzyme then hydrolyzes the chromogenic synthetic substrate. As a result, it produces a yellow coloration. β -glucan also activates β -glucan dependent clotting cascade, however, activation by β -glucan is dependent on the concentration of β -glucan within the sample. In the presence of very high concentration of β -glucan, the lysate will not be activated. The Limulus Color KY Single Test contains a large amount of β -glucan to prevent the activation. This formulation allows β -glucan to be completely inhibited and endotoxin can be assayed specifically⁽¹⁾. The Limulus Color KY Single Test is measured using a Toxinometer[®] to show the activation of product due to color generation when there is a decrease in the transmitted light quantity ratio. In addition, it measures the activation time (Ta), which is the reaction time to reach the preset threshold value. As a result, the endotoxin concentration of the sample is calculated from the relationship between the Ta and endotoxin concentration generated from the standard curve.

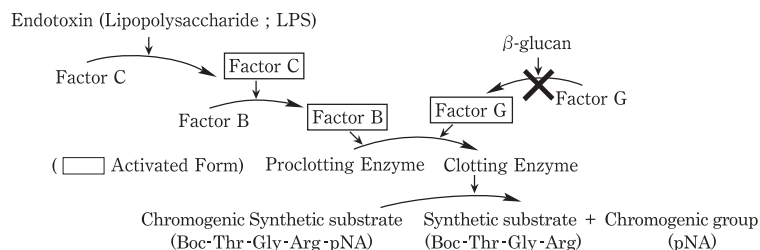
Figure 1 The Color Development Cascade Mechanism of Limulus Color KY Single Test

The × represents that the activation of β -glucan is prevented due to the large amount of β -glucan derivative in the reagent.

[Kit Contents]

1. Limulus Color KY reagent, lyophilized extract of *Limulus polyphemus* amebocytes containing synthetic chromogenic substrate, buffer constituents and (1 → 3)- β -D-glucan derivatives.
..... 25 vials × 0.2 mL

Storage : Store at 2 to 10°C.

2. Control Standard Endotoxin (CSE), lyophilized... 1 vial (500 ng purified LPS) Contains LPS purified from *E. coli* UKT-B including mannitol and glycine as additives.

Storage : Store at 2 to 10°C. The reconstituted CSE solution can be stored at 2-10°C for 1 month. Stored CSE solution must be agitated vigorously with a vortex mixer for more than 1 minute before use.

I. Materials and Equipment not provided

1. Pipettes
2. Test tubes* for dilution of standard solutions
3. Vortex mixer
4. Water for BET*
5. Micropipette tips* (Bio Clean Tip Wako[®] Extend S II Code : 294-35011, Bio Clean Tip Wako[®] 1000 II Code : 298-35031)
6. Toxinometer[®]
7. Aluminum Cap-S* (Code 293-28251)

* : The materials should be endotoxin-free.

NOTE : A test tube and the pipette should be dehydrogenated by 250 °C

and more than 30 min. Water should be Endotoxin-free water. The plastic products such as tips for micropipets should be confirmed that they are no endotoxin contamination and adsorption or interference with the measurement.

II. Preparation of Reagents

1. Preparation of Control Standard Endotoxin (CSE) Solution

- 1) Remove the aluminum cap from the vial.
- 2) Pull up the rubber stopper with tweezers and remove the stopper from the vial.
- 3) Add the appropriate volume of water for BET based on the EU value (EU/vial) printed on the label to make 1000 EU/mL solution.
- 4) Return the rubber stopper on to the vial.
- 5) Agitate the vial vigorously with a Vortex mixer for approximately 2 minutes to dissolve the lyophilized product.

2. Preparation of dilution standards

- 1) Use the CSE solution from the "Preparation of Control Standard Endotoxin (CSE) Solution" to prepare dilution standards.
- 2) Prepare a dilution series, ensure that one dilution does not exceed more than 10-fold.
- 3) Each standard solution should be mixed with a vortex mixer for more than 30 seconds before proceeding to the next dilution in the series.
- 4) The dilution standard solutions should be assayed after preparation.

3. Dilution of samples

Samples often need to be diluted. A diluted sample should be agitated with a vortex mixer for more than 30 seconds before being used in the procedure.

4. pH adjustment of samples

If the pH of a mixture of the sample + Limulus Color KY

Single Test solution is out of the range (between 6.0 and 8.0), adjust the pH range of the sample solution to the appropriate range. This can be achieved with a diluted sodium hydroxide or diluted hydrochloric acid solution.

III. Determination using Toxinometer[®] Measurement System^{13), 14), 15), 16)}

1. Preparation of Reagents

Prepare reagents as stated in the "Preparation of Reagents" section.

2. Toxinometer[®] Parameters

Set the Toxinometer[®] parameters as shown :

Target temperature : 37°C
Threshold value (Th) : 94.9%
Count : 3
Wait time^{*1} : 5 minutes
Measurement time^{*2} : 60 minutes

^{*1} : When measuring at high concentrations, it is recommended to shorten the Wait time and change the setting to "Use the first count as Reaction time."

^{*2} : Please note that the higher the sensitivity, the more time may be required for the endotoxin to be measured by the Toxinometer[®].

3. Assay procedure

- 1) Remove the appropriate number of Limulus Color KY Single Test vials out of the kit box and place each in the test tube rack
- 2) Remove the cap by using the cap opener provided in the kit box
- 3) If the lyophilized product is attached to the rubber stopper, gently tap the bottom of the vial on a firm surface before opening the vial. Small amounts of lyophilized product attached to the rubber stopper will not affect testing.
- 4) Use cap opener to pull up rubber stopper releasing the vacuum. If

necessary, remove the rubber stopper and cover vial with an endotoxin-free aluminum cap.

- 5) Dispense 0.2 mL of water for BET (negative control), standard solutions, and/or samples into each of the lyophilized Limulus Color KY Single Test vials.
- 6) Return the rubber stopper or the aluminum cap on to the single test vial.
- 7) Mix the single test vial with a vortex mixer for 5 seconds to dissolve. Check to dissolved.
- 8) Place each vial into the Toxinometer[®] to be measured by the Toximaster[®] software

4. Standard curve and Data analysis

The endotoxin concentration of the sample is based on the sample's activation time using the generated standard curve.

A standard curve will be generated by taking the logarithm of the endotoxin concentration on the x axis and the logarithm (log-log) or twice the logarithm (log-loglog) of the activation time on the y axis. Linear or quadratic regression can be used for regression calculation. If the range of the standard curve is wide, improvement of the standard curve fitting can be expected by using twice the logarithm of activation time or quadratic regression.

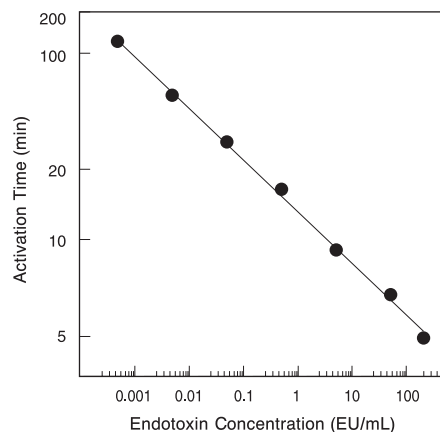
Figure 2 : An example of standard curve of endotoxin measured with Limulus Color KY Single Test using a Toxinometer[®].

This figure shows that a high correlation standard curve can be obtained over a wide concentration range.

[Warnings and General Precautions]

- Do not use this reagent for any purposes other than BET.
- This product is for research use only, not for diagnostic use.
- Exercise caution when handling the reagent, because toxicity of LAL is unknown.
- Even a trace amount of the Limulus Color KY Single Test reagent, which contains a large amount of β -glucan derivative, may activate regular LAL reagent that is sensitive to β -glucan. Be careful not to contaminate regular LAL with the Limulus Color KY Single Test reagent.
- Be careful not to cause endotoxin contamination and always use aseptic techniques.
- Discard product if the color has changed to yellow or any insoluble matter is found after dissolution.

• Do not use reagents past their expiration dates.

[References]

- 1) Bang, F.B. *Bull Johns Hopkins Hosp.* **98**, 325-351 (1956).
- 2) Levin, J. and Bang, F.B. *Bull Johns Hopkins Hosp.* **115**, 265-274 (1964).
- 3) *The United States Pharmacopeia 41st, the National Formulary 36*, 6011, U.S. Pharmacopeial Convention Inc. (2018).
- 4) *The Japanese Pharmacopoeia, Seventeenth Edition Handbook*, B-473-B-492 (2016).
- 5) Kakinuma, J., Asano, T., Torii, H., and Sugino, Y. : *Biochem. Biophys. Res. Commun.*, **101**, 434-439 (1981).
- 6) Nakamura T, Morita T, Hiranaga M, Miyata T, Iwanaga S. : *Nihon Saikingaku Zasshi*, **38**, 781-803 (1983).
- 7) Iwanaga, S., Morita, T., Miyata, T., Nakamura, T., Hiranaga, M. and Ohtsubo, S. : *Bacterial Endotoxin : Chemical, Biological and Clinical Aspects*, edited by Homma, J. Y., Kanegasaki, S., Luderitz, O., Shiba, T. and Westphal, O., 365-382, Verlag Chemie (1984).
- 8) Peason, F. C., Weary, M. and Bohon, J. : *Endotoxins and Their Detection with the Limulus Amebocyte Lysate Test*, edited by Watson, S. W., Levin, J. and Novitsky, T. J., 247-260, Alan R. Liss, Inc. (1982).
- 9) Westphal, O., Luderiz, O., Eichenberger, E. and Keiderling, W. : *Z. Naturforsch.*, **7b**, 536-548 (1952).
- 10) Akama, K., Kuratsuka, K. and Homma, R. Kanoh, S., Niwa, M., Iwanaga, S. and Nakahara, C. : *Bacterial Endotoxin : Chemical, Biological and Clinical Aspects*, edited by Homma, J. Y., Kanegasaki, S., Luderitz, O., Shiba, T. and Westphal, O., 395-407, Verlag Chemie (1984).
- 11) Tsuchiya M, Takaoka A, Tokioka N, Matsuura S. *Nihon Saikingaku Zasshi*. **45**, 903-911 (1990).
- 12) Oishi, H., Hatayama, Y., Shiraishi, H., Yanagisawa, K., Sakata, Y. : *Yakugakuzasshi*, **105**, 300 (1985).
- 13) Oishi, H., Takaoka, A., Hatayama, Y., Matsuo, T. and Sakata, Y. : *J. Parenter. Sci. Technol.*, **39**, 194-199 (1985).
- 14) Oishi, H., Fusamoto, M., Hatayama, Y., Tsuchiya, M., Takaoka, A. and Sakata, Y. : *Chem. Pharm. Bull.*, **36**, 3012-3019 (1988).
- 15) Oishi, H., Wada S., Shiraishi H., Koshindo T., Tsujino R. : The 116th Annual Meeting the Pharmaceutical Society of Japan, 100, (1996).

[Storage] 2 ~ 10°C

[Expiration date] Indicated on the label

[Package] 25 tests

Manufactured by
FUJIFILM Healthcare Manufacturing Corporation

2613-2, Ogohara, Komono-Cho, Mie-Gun, Mie 510-1222, Japan

Distributed by
FUJIFILM Wako Pure Chemical Corporation

1-2, Doshomachi 3-Chome, Chuo-Ku, Osaka 540-8605, Japan
Telephone : + 81-6-6203-3741
Facsimile : + 81-6-6201-5964
<http://ffwk.fujifilm.co.jp>

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