LBIS™ Human IFN- γ ELISA Kit

Please, read this instruction carefully before use.

1. Intended use

Interferon-gamma (IFN- γ) is the only type II interferon consisting of 143 amino acid residues (human), which is structurally different from type I interferons such as IFN- α/β . IFN- γ is secreted predominantly by T cells and NK cells and functions as a homodimer. It is known to have various biological activities, including antiviral effect, cell growth inhibition, anti-tumor effect, alteration of macrophage activity, promotion of NK cell activity, immunomodulatory effects, modulation of the immune response, modulation of the differentiation process, and is attracting attention in various research fields.

This kit is a sandwich ELISA system for quantitative measurement of human IFN- γ . It can measure IFN- γ in human serum (plasma) specifically and with high sensitivity. This product is intended for research use only. Not for use in diagnostics procedures.

2. Assay Principle

In FUJIFILM Wako Pure Chemical Corporation's LBIS™ Human IFN- γ ELISA Kit, biotinylated anti-IFN- γ antibody and standards or samples are incubated in monoclonal antibody-coated wells to capture IFN- γ . After 90 minutes incubation and washing, HRP (horse radish peroxidase)-conjugated streptavidin is added and incubated further for 30 minutes to bind with captured IFN- γ . After washing, bound HRP-conjugated streptavidin is reacted with a TMB Solution for 20 minutes, and reaction is stopped by addition of acidic solution, and absorbance of yellow product is measured spectrophotometrically at 450 nm. The absorbance is nearly proportional to IFN- γ concentration. The standard curve is prepared by plotting absorbance against standard concentrations. IFN- γ concentrations in unknown samples are determined using this standard curve.

3. Performance Characteristics

- Precision of assay (Within assay variation)
2 serum samples, 5 replicates assay

n/ID	Sample1 (pg/mL)	Sample2 (pg/mL)
1	26.1	1.41
2	26.4	1.42
3	26.2	1.44
4	26.3	1.46
5	26.4	1.43
mean	26.3	1.43
SD	0.13	0.019
CV (%)	0.50	1.3

- Reproducibility (Between assay variation)
3 serum samples, 4 days, 3 replicates assay

Day/ID	Sample1 (pg/mL)	Sample2 (pg/mL)	Sample3 (pg/mL)
0	29.1	4.90	0.742
1	29.9	4.69	0.723
2	31.2	4.87	0.785
3	29.5	4.99	0.772
mean	29.9	4.86	0.756
SD	0.92	0.13	0.028
CV (%)	3.1	2.6	3.7

- Recovery test (Standard human IFN- γ was added in 3 concentrations to serum or plasma sample and assayed.) The recoveries were 86.8% - 108%


Serum sample

Added (pg/mL)	Found (pg/mL)	Recovered (pg/mL)	Recovery (%)
–	0.516	–	–
1.00	1.60	1.08	108
4.00	4.50	3.98	99.5
16.0	16.4	15.9	99.4

Plasma sample (heparin)

Added (pg/mL)	Found (pg/mL)	Recovered (pg/mL)	Recovery (%)
–	0.618	–	–
1.00	1.62	1.00	100
4.00	4.20	3.58	89.6
16.0	14.5	13.9	86.8

- Dilution test (Serum sample and plasma sample were serially diluted by 3 steps.)

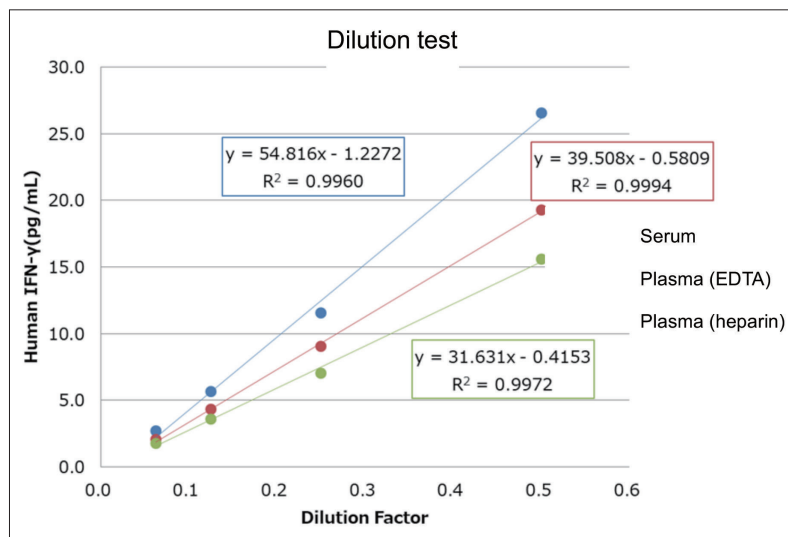4. Technical Tips/Precautions

- ELISA can be affected by assay environment. Ensure that room temperature at places for assay operation and incubation is strictly controlled at 20°C - 25°C. Avoid performing assay under the airstream velocity over 0.4 m/sec. and the humidity less than 30%, seal the well plate with a plate seal and place the well plate in an incubator or a styrofoam box in each step of incubation.
- To avoid denaturation of the coated antibody, do not let the plate go dry.
- In order to avoid dryness of wells, contamination of foreign substances and evaporation of dispensed reagents, never forget to cover the well plate with a plate seal supplied, during incubation.
- Be careful to avoid any contamination of assay samples and reagents. We recommend the use of disposal pipette tips, and 1 tip for 1 well.
- As the Antibody-coated Plate is module type of 8 wells $\times$ 12 strips, each strip can be separated by cutting the cover sheet with a knife and used independently.
- The TMB Solution should be almost clear pale yellow before use. It turns blue during reaction and gives yellowish color after addition of Stop Solution. Greenish color means incomplete mixing.
- Time the reaction from the pipetting of the reagent to the first well.
- Prepare a standard curve for each assay.
- For professional use only, beginners are advised to use this kit under the guidance of experienced person. In manual operation, proficiency in pipetting technique is recommended.
- Wear protective gloves, clothing, eye, and face protection.
- Avoid contact of skin and mucous membranes with kit reagents or specimens. If any reagents come in contact with eyes, skin, or mucous membranes, wash with copious amounts of water and contact a physician.
- Do not drink, eat, or smoke in a place where this kit is used.
- The materials must not be pipetted by mouth.
- Do not use reagents with different lot numbers together.
- The reagents are prepared to give accurate results only when used in combination within the same box. Therefore, do not combine the reagents from kits with different lot numbers. Even if the lot number is the same, it is best not to mix the reagents with those that have been preserved for some period.
- Handle the sample with proper care, being aware that the sample may have an infection risk. This kit contains animal-derived ingredients.
- Residual samples and used tips should be sterilized before disposal.
- Dispose consumable materials and unused contents in accordance with applicable regional/national regulatory requirements.

5. Reagents supplied

Components	Use Status	Amount
(A) Antibody-coated Plate	Use after washing.	96 wells/1 plate
(B) Human IFN- γ Standard	Freeze-dried. Use after reconstitution.	1 bottle
(C) Buffer	Ready for use.	60 mL/1 bottle
(D) Biotin-conjugated Antibody	Freeze-dried. Use after reconstitution.	1 bottle
(E) Peroxidase-conjugated Streptavidin Solution	Concentrated. Use after dilution.	100 μ L/1 bottle
(F) TMB Solution	Ready for use.	12 mL/1 bottle
(H) Stop Solution	Ready for use.	12 mL/1 bottle
(I) Wash Solution (10 $\times$)	Concentrated. Use after dilution.	100 mL/1 bottle
(J) Plate Seal	—	4 sheets

[Storage and Stability]

[(A) Antibody-coated Plate]

If seal is not removed, put the strip back in a plastic bag with zip lock originally used for well-plate container and store at 2°C - 10°C.

[(B) Human IFN- γ Standard]

The reconstituted standard solution (original standard solution) should be used within 2 weeks if stored in a refrigerator. Dispose the remaining diluted standard solutions after use.

[(C) Buffer] & [(F) TMB Solution]

Use only volume you need for your assay. Remaining reagents should be stored at 2°C - 10°C closing the cap tightly. Once opened, we recommend using as soon as possible to avoid influence by environmental condition.

[(D) Biotin-conjugated Antibody]

The reconstituted (D) solution should be used within 2 weeks if stored in a refrigerator. Dispose the remaining diluted working solution.

[(E) Peroxidase-conjugated Streptavidin Solution]

Remaining working solution (already diluted) should be disposed. The rest of the undiluted solution (unused) : close the cap tightly and store at 2°C - 10°C. Once opened, we recommend using as soon as possible to avoid influence by environmental condition.

[(H) Stop Solution]

Close the cap tightly and store at 2°C - 10°C.

[(I) Wash Solution (10 $\times$)]

The rest of undiluted solution (unused) : close the cap tightly and store at 2°C - 10°C. Dispose any remaining diluted buffer.

6. Equipments or supplies required but not supplied ☐ Use as a check box

☐ Deionized water (or Distilled water) ☐ Test tubes for preparation of standard solution series. ☐ Glassware for dilution of Wash Solution (10 $\times$) (a graduated cylinder, a bottle) ☐ Pipettes (disposable tip type). One should be able to deliver 50 μ L precisely and another for 100 μ L - 1000 μ L. ☐ Syringe-type repeating dispenser like Eppendorf multipette plus which can dispense 100 μ L. ☐ Paper towel to remove washing buffer remaining in wells. ☐ A vortex-type mixer. ☐ A shaker for 96 well-plate (500 rpm - 1200 rpm) ☐ An automatic washer for 96 well-plate (if available), or a wash bottle with a jet nozzle. ☐ A 96 well-plate reader (450 nm $\pm$ 10 nm, 620 nm : 600 nm - 650 nm) ☐ Software for data analysis.

7. Preparation of Samples

- This kit is intended to measure IFN- γ in human serum or plasma.
- We recommend using EDTA-2Na, EDTA-2K or heparin as anticoagulant.
- Samples should be immediately assayed or stored below -35°C until assay. Before starting assay, shake thawed samples sufficiently. Do not repeat freeze-and-thaw cycles.
- Centrifuge the sample to remove turbidity or insoluble matters where necessary before use for the assay
- If presence of interfering substance is suspected, examine by dilution test at more than 2 points.
- Dilution of a sample should be made in a PP or PE test tube using Buffer prior to adding them to wells

8. Preparation of Reagents

- ◆ Bring all reagents of the kit to room temperature (20°C - 25°C) before use.
- ◆ Prepare reagent solutions in appropriate volume for your assay. Do not store the diluted reagents.

[Concentrated reagents]

[(B) Human IFN- γ Standard]

Reconstitute the (B) Human IFN- γ standard with purified water in the amount specified in "Reconstitution of standard"* to prepare a original standard solution (750 pg/mL).


Then prepare using the kit-attached (C) buffer that has been allowed to warm up to room temperature.

*Please confirm "Reconstitution of standard" from this product page. Since the amount of purified water added varies by lot, please be sure to check each lot.

Below is an example of preparing each standard solution.

Volume of standard solution	(C) Buffer	Concentration
Original standard solution 50 μ L	450 μ L	75.0 pg/mL
75.0 pg/mL solution 200 μ L	300 μ L	30.0 pg/mL
30.0 pg/mL solution 200 μ L	300 μ L	12.0 pg/mL
12.0 pg/mL solution 200 μ L	300 μ L	4.80 pg/mL
4.80 pg/mL solution 200 μ L	300 μ L	1.92 pg/mL
1.92 pg/mL solution 200 μ L	300 μ L	0.768 pg/mL
Blank	300 μ L	0

Standard Curve

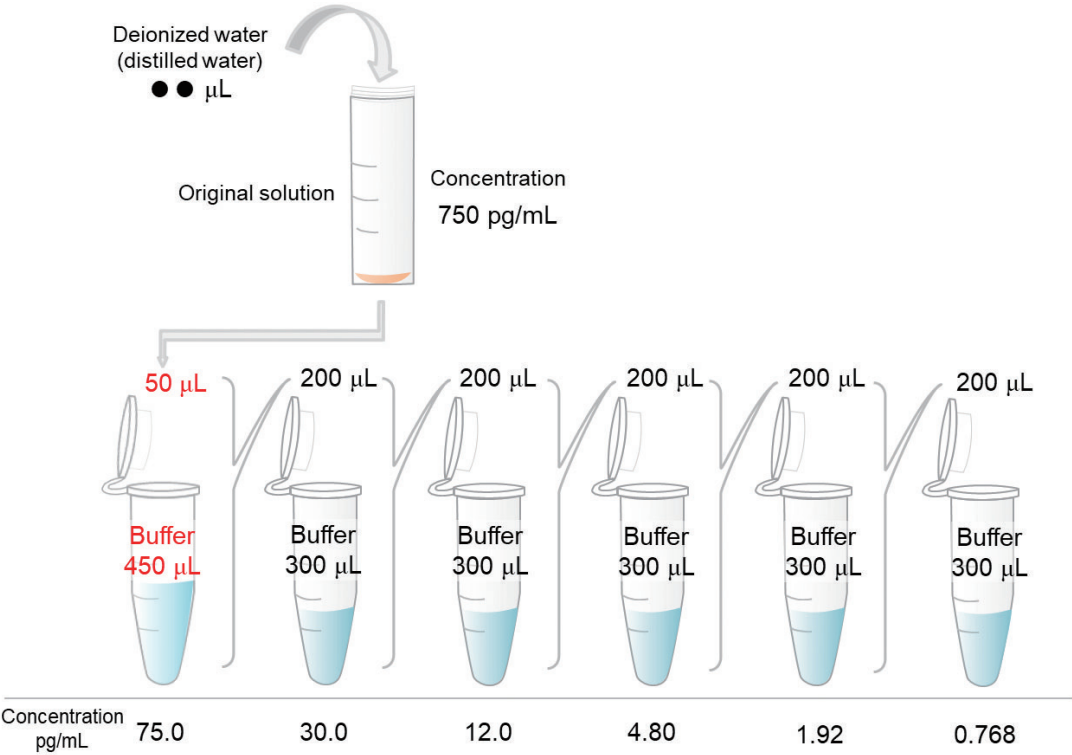[(D) Biotin-conjugated Antibody]

Add 100 μ L of deionized water (or distilled water) of room temperature to (D) Biotin-conjugated Antibody (Freeze-dried) (Original (D) solution). Prepare working solution by dilution of original (D) solution with the (C) Buffer to **1 : 100**. 10 mL of the diluted solution is enough for 96 wells.

[(E) Peroxidase-conjugated Streptavidin Solution]

Prepare working solution by dilution of (E) with the (C) Buffer to **1 : 100**. 10 mL of the diluted solution is enough for 96 wells.

[(I) Wash Solution (10 $\times$)]

Dilute 1 volume of the concentrated (I) Wash Solution (10 $\times$) to **10 volumes** with deionized water (or distilled water) to prepare working solution. Example : 100 mL of concentrated (I) Wash Solution (10 $\times$) and 900 mL of deionized water (or distilled water).


9. Assay Procedure

Remove the cover sheet of the Antibody-coated Plate after bringing up to room temperature.

- (1) Wash the Antibody-coated Plate (A) by filling the wells with washing buffer and discard 4 times (*①), then strike the plate upside-down onto folded several sheets of paper towel to remove residual buffer in the wells. (*⑤)
 - (2) Pipette 50 μ L of Biotin-conjugated Antibody solution to all wells, and shake the plate gently on a plate shaker (*②).
 - (3) Pipette 100 μ L of standard solution to the wells designated for standards.
 - (4) Pipette 50 μ L of (C) Buffer the sample wells, and then add 50 μ L of each sample to the wells.
 - (5) Stick a Plate Seal (*③) on the plate and shake the plate gently on a plate shaker (500 rpm).
Incubate for 90 minutes at 20°C - 25°C.
 - (6) Discard the reaction mixture and rinse wells as step (1).
 - (7) Pipette 100 μ L of Peroxidase-conjugated Streptavidin Solution to all wells and shake as step (4).
 - (8) Stick a Plate Seal (*③) on the plate and incubate the plate for 30 minutes at 20°C - 25°C.
 - (9) Discard the reaction mixture and rinse wells as step (1).
 - (10) Pipette 100 μ L of TMB Solution to wells and shake as step (4).
 - (11) Stick a Plate Seal (*③) on the plate and incubate the plate for 20 minutes at 20°C - 25°C.
 - (12) Add 100 μ L of the Stop Solution to all wells and shake as step (4).
 - (13) Measure the absorbance of each well at 450 nm (reference wavelength, 620 nm (*④)) immediately using a plate reader.
- *Refer to 13. Summary of Assay Procedure for notes on *①, *②, *③, *④, and *⑤.

10. Calculations

- (1) Create a standard curve by reducing the data using computer software capable of generating 3rd order regression curve, or 4 or 5 parameters. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the human IFN- γ concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.
- (2) Using the standard curve, read the IFN- γ concentration of samples at its absorbance, and multiply the assay value by dilution factor if the sample has been diluted. The standard operation method is 2-fold dilution.

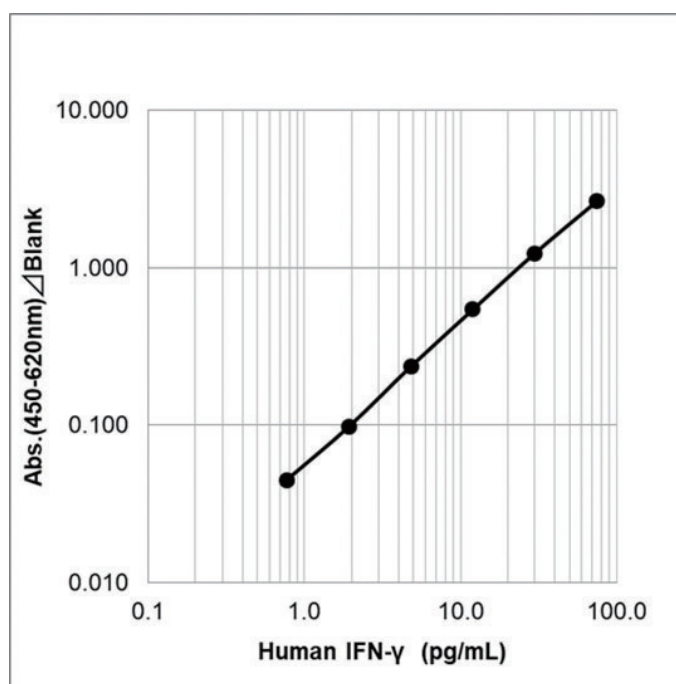11. Calibration

The measured value obtained by this kit can be converted to the unit concentration based on NIBSC/WHO standard IFN- γ (Code : 82/587) by multiplying by the conversion factor specified in a separate sheet according to the following calculation formula.

NIBSC/WHO (82/587) unit concentration (U/mL) = conversion factor (*1) $\times$ measured IFN- γ value (pg/mL) (*2)

(*1) The exchange rate is different depending on the lot. Please see the specified conversion factor in the annex document.

(*2) Use the assay value that multiplied by dilution rate.

12. Trouble Shooting

- Low absorbance in all wells

Possible explanations :

- 1) The standard or samples might not be added.
- 2) Reagents necessary for coloration such as Biotin-conjugated Antibody, Peroxidase-conjugated Streptavidin Solution, or TMB Solution might not be added.
- 3) Wrong reagents related to coloration might have been added. Wrong dilution of Biotin-conjugated Antibody or Peroxidase-conjugated Streptavidin Solution.
- 4) Contamination of enzyme inhibitor(s).
- 5) Influence of the temperature under which the kits had been stored.
- 6) Excessive hard washing of the well plate.
- 7) Addition of chromogenic substrate reagent soon after taking out from a refrigerator might cause poor coloration owing to low temperature.

- Blank OD was higher than that of the lowest standard concentration (0.768 pg/mL).

Possible explanations :

Improper or inadequate washing. (Change washing frequency from 4 times to 5 - 8 times at the constant stroke after the reaction with Peroxidase-conjugated Streptavidin Solution.)

- High coefficient of variation (CV)

Possible explanation :

- 1) Improper or inadequate washing.
- 2) Improper mixing of standard or samples.
- 3) Pipetting at irregular intervals.

- Q-1 : Can I divide the plate to use it for the other testing?

A-1 : Yes, cut off the clear seal on the plate with cutter along strip. Put the residual plate, which is still the seal on, in a refrigerator soon.

- Q-2 : I found 96 well-plate is empty when I opened the box.

A-2 : As this kit is dried type, not preservation stabilizer is added.

13. Summary of Assay Procedure ☐ : Use as a check box

*First, read this instruction manual carefully and start your assay after confirmation of details.

☐ Bring the well-plate and all reagents back to **20°C - 25°C for 2 hours**.

☐ (I) Wash Solution (10×) must be diluted to **10 times** by deionized water (or distilled water) that returned to 20°C - 25°C.

☐ Standard solution dilution example

Reconstitute the (B) Human IFN- γ standard with purified water in the amount specified in "Reconstitution of standard"* to prepare a original standard solution (750 pg/mL).

Then prepare using the kit-attached (C) buffer that has been allowed to warm up to room temperature.

*Please confirm "Reconstitution of standard" from this product page. Since the amount of purified water added varies by lot, please be sure to check each lot.

Preparation of diluted IFN- γ standard solutions :


Conc. (pg/mL)	75.0	30.0	12.0	4.80	1.92	0.768	0
Std. sol. (μ L)	50	200*	200*	200*	200*	200*	—
Buffer (μ L)	450	300	300	300	300	300	300

*One rank higher standard.

☐ Add deionized water (or distilled water) of room temperature to (D) Biotin-conjugated Antibody (Freeze-dried) (Original (D) solution). Prepare working solution by dilution of original (D) solution to **1 : 100** with the (C) Buffer returned to 20°C - 25°C.

- ☐ Antibody-coated Plate
 - ☐ ↓ Washing 4 times (*①), (*⑤)
 - ☐ Biotin-conjugated Antibody solution (Diluted) 50 μ L
 - ☐ ↓ Shaking (*②)
 - ☐ (Diluted) samples, or Standards 100 μ L
 - ☐ ↓ Shaking (500 rpm), Incubation for 90 minutes at 20°C - 25°C.
 - ☐ Dilute Peroxidase-conjugated Streptavidin to **1 : 100** with the (C) Buffer returned to 20°C - 25°C.
(This should be prepared during incubation.)
 - ☐ ↓ Washing 4 times (*①), (*⑤)
 - ☐ Peroxidase-conjugated Streptavidin Solution (Diluted) 100 μ L
 - ☐ ↓ Shaking (*②), Incubation for 30 minutes at 20°C - 25°C. (Standing (*③))
 - ☐ ↓ Washing 4 times (*①), (*⑤)
 - ☐ TMB Solution 100 μ L
After dispensing, the color turns to blue depending on the concentration.
 - ☐ ↓ Shaking (*②), Incubation for 20 minutes at 20°C - 25°C. (Standing (*③))
 - ☐ Stop Solution 100 μ L
After dispense, the color turns to yellow depending on the concentration.
 - ☐ ↓ Shaking (*②), Immediately shake.
 - ☐ Measurement of absorbance (450 nm, Ref 620 nm (*④)) immediately.
Ref. wave cancels the dirt in the back of plate.
- *① After dispensing wash buffer to wells, lightly shake the plate on your palm for 10 sec and remove the buffer. Guideline of washing volume : 300 μ L/well for an automatic washer and for a pipette if the washing buffer is added by pipette. In case of washing by using 8 channel pipettes, sometimes the background tends to be high. If so, change washing frequency from 4 times to 5 - 8 times at the constant stroke after the reaction with Peroxidase-conjugated Streptavidin Solution.
Standard of plate-washing pressure : 5 mL/min - 25 mL/min. (Adjust it depending on the nozzle's diameter.)
- *② Guideline of shaking : 500 rpm - 1200 rpm 3 times for 10 seconds.
- *③ Seal the plate during the reaction after shaking. Peel off the protective paper from the seal and stick the seal on the plate. Do not reuse the plate seal used once.
- *④ 600 nm - 650 nm can be used as reference wavelength.
- *⑤ After removal of wash buffer, immediately dispense the next reagent.

Assay Worksheet

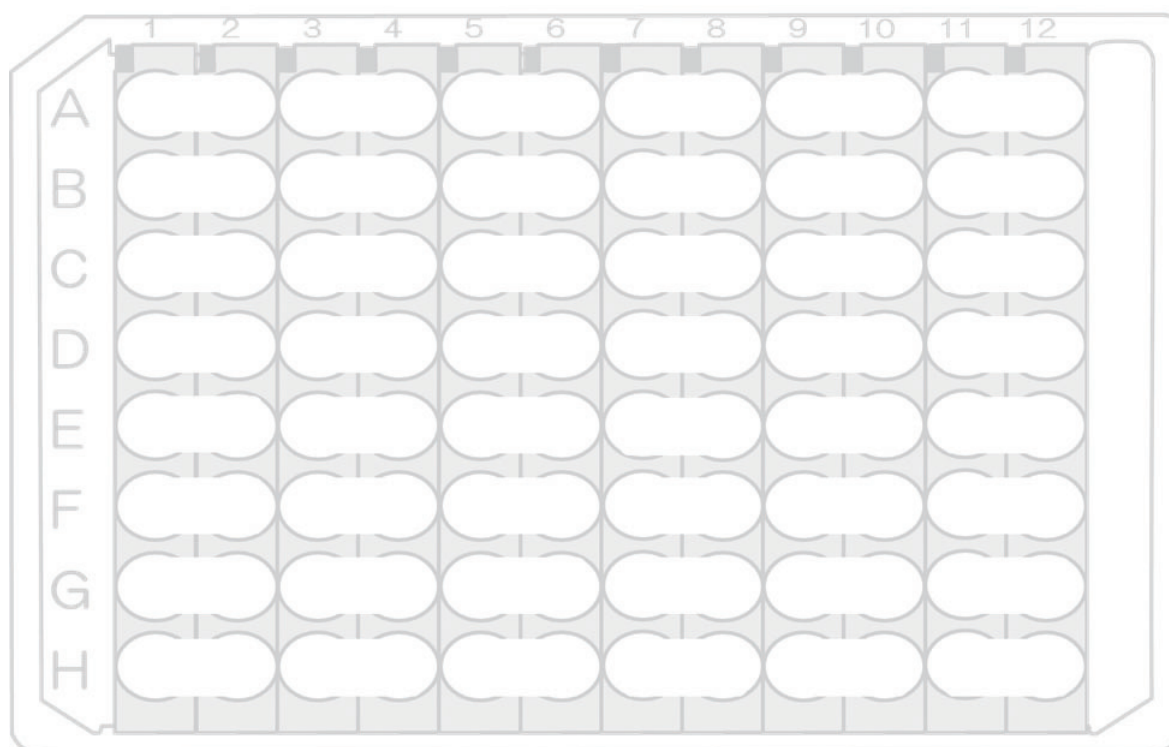14. Storage and Expiration

The complete kit is stored at 2°C - 10°C (Do not freeze). Opened reagents should be used as soon as possible to avoid less than optimal assay performance caused by storage environment.

LBIS™ Human IFN- γ ELISA Kit

[Storage] Store the kit at 2°C - 10°C (Do not freeze)

[Expiration date] Indicated on the container

[Package] For 96 tests

[Cat #] 298-87301

FUJIFILM Wako Pure Chemical Corporation

1-2, Doshomachi 3-Chome, Chuo-Ku, Osaka 540-8605, Japan
Telephone : + 81-6-6203-3741
Facsimile : + 81-6-6201-5964
<http://flwk.fujifilm.co.jp>

FUJIFILM Wako Chemicals U.S.A. Corporation

1600 Bellwood Road
Richmond, VA 23237
U.S.A.
Telephone : + 1-804-271-7677
Facsimile : + 1-804-271-7791
<http://www.wakousa.com>

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