




Solvents for Liquid Chromatography

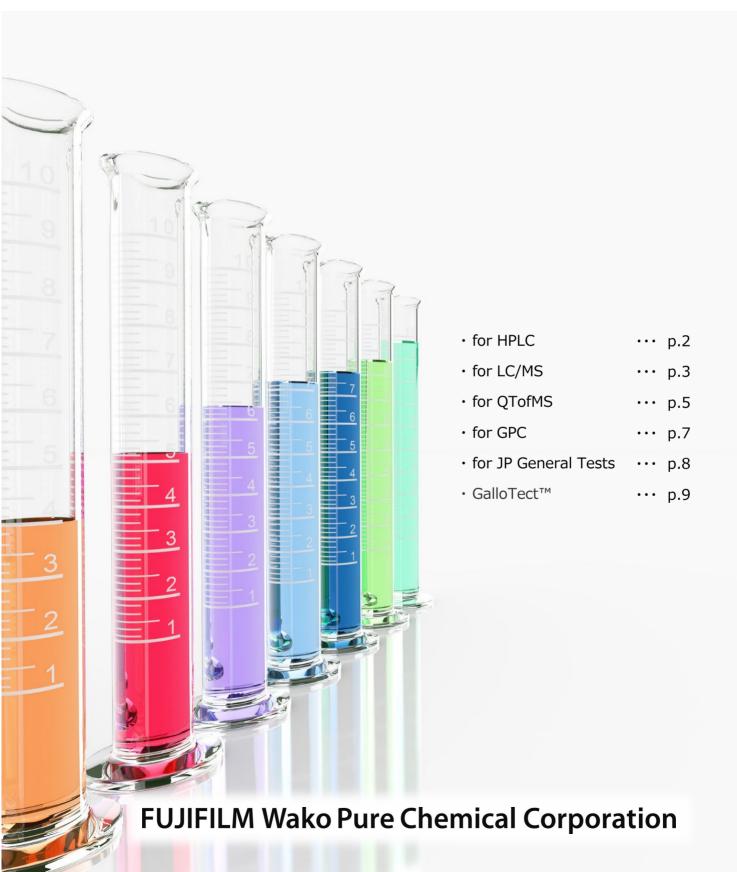


LC solvents of FUJIFILM Wako Pure Chemical Corporation

As various chromatograph detectors become more improved, the quality demanded in the mobile phase (solvent) has become higher and more diversified. We have a lineup of solvents for high-performance liquid chromatograph (HPLC) that can be used in LC analysis in general, solvents for high-purity LC/MS, and solvents for QTofMS, which is suitable for ultrahigh sensitivity MS analysis. We also offer optimal solvents for use in GPC, and solvents that can be used for Japanese Pharmacopoeia tests.

List of testing requirement for LC solvents

The proof test items vary even for the same solvent, depending on the detector in which it is expected to be used or the application.

The solvents under each standard have-passed the testing requirement that are suited to the "expected applications," in addition to the quality test as chemical substances.

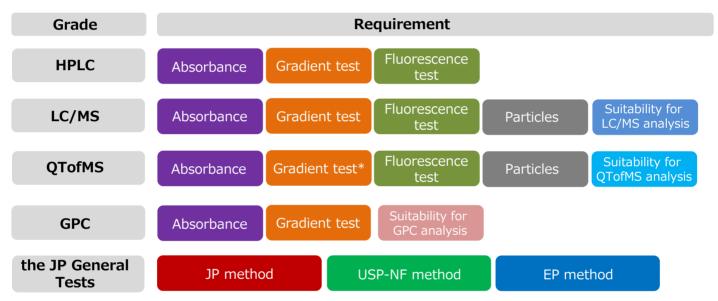

^{*}The gradient test for QTofMS was conducted as a QTofMS-compatible test.

Description of testing requirement

A test to confirm that the peak hight that interferes with the analysis is lower than **Absorbance** a certain standard level at a specific wavelength. A test to confirm that the fluctuation of the baseline is lower than at the standard Gradient test level when it is used in gradient elution. A test that compares the quantity of impurities with fluorescence that cause the background Fluorescence noise or deterioration in detection sensitivity when it is used in fluorescence measurement test with the reference sample by the fluorescence intensity. **Particles** A test to confirm the number of particles that may affect the column life. A test to confirm that the background noise is suppressed to a low level when it is used in Suitability for analysis using a mass spectrometer. LC/MS analysis A test to confirm that that baseline fluctuation is lower than the standard level by connecting Suitability for a column to UHPLC and conducting gradient measurement using QTofMS. QTofMS analysis A test to confirm the absorbance after heating with consideration of solvent discoloration

GPC analysis

A test to commit the absorbance after heating with consideration of solvent discoloration and analysis.

for High Performance Liquid Chromatography (HPLC)

Absorbance

Gradient test

Fluorescence test

We offer various solvents of quality appropriate for use in high-performance liquid chromatograph.

Product list (Solvents for HPLC)

| Product Number | Product Name | Package Size |
|----------------|---|--------------|
| 014-08681 | Acatona | 1 L |
| 010-08683 | Acetone | 3 L |
| 019-08631 | Acatonituila | 1 L |
| 015-08633 | Acetonitrile | 3 L |
| 025-06691 | Benzene | 1 L |
| 023-10801 | 1-Butanol | 1 L |
| 024-12771 | t-Butyl Methyl Ether | 1 L |
| 033-08631 | Chlavafavva | 1 L |
| 039-08633 | Chloroform | 3 L |
| 031-20531 | Chlarafarma Amariana addad | 1 L |
| 037-20533 | Chloroform, Amylene added | 3 L |
| 033-08511 | Cyclebayana | 1 L |
| 039-08513 | Cyclohexane | 3 L |
| 042-16691 | 1,4-Dioxane | 1 L |
| 046-18671 | o-Dichlorobenzene | 1 L |
| 042-20621 | N. N. Dimothylformamida | 1 L |
| 048-20623 | <i>N,N</i> -Dimethylformamide | 3 L |
| 046-16971 | Distilled Water | 1 L |
| 042-16973 | Distilled Water | 3 L |
| 136-06751 | Dichloromethane | 1 L |
| 132-06753 | Dictionormediane | 3 L |
| 056-03341 | [Floring (00 F) | 1 L |
| 052-03343 | Ethanol (99.5) | 3 L |
| 057-03371 | Ethyl Acetate | 1 L |
| 053-03373 | Ltityi Acetate | 3 L |
| 085-06991 | 1,1,1,3,3,3-Hexafluoro-2-propanol | 100 mL |
| 087-06995 | 1,1,1,3,3,3-1 lexalition 0-2-pi oparior | 500 mL |
| 084-03421 | Hexane | 1 L |
| 080-03423 | I lexalle | 3 L |
| 085-03691 | Heptane | 1 L |
| 138-06473 | Methanol | 1 L |
| 132-06471 | r-ica ianoi | 3 L |
| 162-13461 | 1-Propanol | 1 L |
| 165-09161 | 2-Propanol | 1 L |
| 161-09163 | 2 1 1 Ο Ράτιοι | 3 L |
| 200-19391 | Tetrahydrofuran, with Stabilizer | 1 L |
| 206-19393 | redanyarolaran, with Stabilizer | 3 L |
| 209-06811 | Tetrahydrofuran, Stabilizer Free | 1 L |
| 205-06813 | redanyarolaran, Stabilizer Free | 3 L |
| 207-06731 | 2,2,4-Trimethylpentane | 1 L |
| 209-06791 | Toluene | 1 L |

for LC/MS

Absorbance

Gradient test

Fluorescence

Particles

Suitability for LC/MS analysis


These are dedicated solvents for LC/MS for which LC/MS compatibility has been confirmed. They are also recommended for ultrahigh-performance liquid chromatograph (UHPLC), as they are guaranteed for particles.

A glass bottle with special treatment to suppress the elution of impurities is used as the container*.

*Patent No. WO/2005/047884

Point for high recommendation 1

The cap material for LC/MS solvents was changed from aluminum to a special resin.

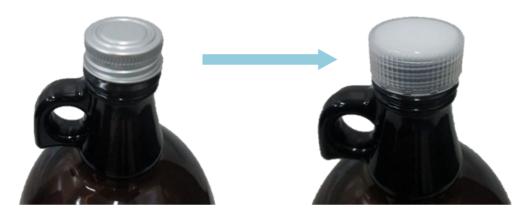


While the conventional aluminum caps posed a risk of aluminum dust adhering to the bottle mouth and contaminating the solvent, the new resin cap has no dust generation or component elution from the material. It maintains the high-quality state of the solvent.

Point for high recommendation 2

The quality was improved further by a review on the manufacturing method. The baseline is suppressed at an even lower level compared to the conventional products.

Comparison of baselines for acetonitrile (comparison of our products)

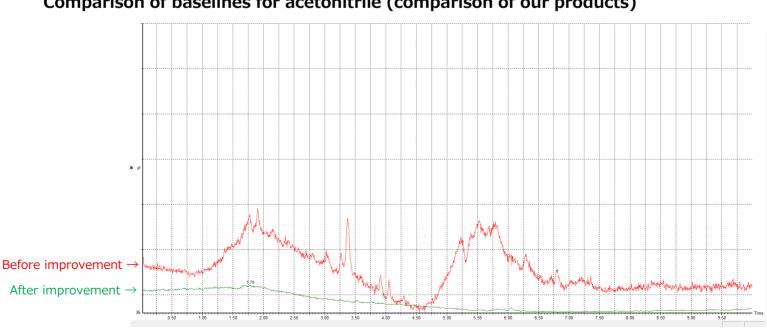

Example of product specifications (acetonitrile)

| | DECLIDENCAL | CDECIFICATION |
|---|--|--------------------|
| | REQUIREMENT | SPECIFICATION |
| Assay (%) | | min. 99.9 |
| Density (20℃ | C) [g/mL] | $0.780{\sim}0.783$ |
| Refractive ind | $ex n_D^{20}$ | $1.343 \sim 1.346$ |
| Water [%] | | max. 0.03 |
| Residue after | evaporation [%] | max. 0.001 |
| Acidity (as Cl | H₃COOH) [%] | max. 0.001 |
| Ammonium (| NH ₄) [ppm] | max. 0.3 |
| Peroxides (as | GH ₂ O ₂) [ppm] | max. 5 |
| Substances reducing permanganate (as O) | | to pass test |
| Gradient test | | to pass test |
| Particles (0.5 [number of p | µm or more) particles per mL] | max. 100 |
| Absorbance | 200nm | max. 0.04以下 |
| | 210nm | max. 0.03以下 |
| | 220nm | max. 0.02以下 |
| | 230nm | max. 0.01以下 |
| | 240nm | max. 0.005以下 |
| Fluorescence | test | to pass test |
| Suitability for LC/MS analysis | | to pass test |

Product list (Solvents for LC/MS)

| Product Number | Product Name | Package Size |
|----------------|----------------------------------|--------------|
| 012-19851 | Acetonitrile | 1 L |
| 018-19853 | Acetonitrile | 3 L |
| 050-09221 | Ethanal (00 E) | 1 L |
| 056-09223 | Ethanol (99.5) | 3 L |
| 062-04721 | | 1 L |
| 068-04723 | 0.1vol% Formic Acid-Acetonitrile | 3 L |
| 138-14521 | Mothanal | 1 L |
| 134-14523 | Methanol | 3 L |
| 168-25531 | 2 Dyamanal | 1 L |
| 164-25533 | 2-Propanol | 3 L |
| 214-01301 | Liltunnium Water | 1 L |
| 210-01303 | Ultrapure Water | 3 L |

Product list (Mobile phase additives for LC/MS)

| Product Number | Product Name | Package Size |
|----------------|------------------------|--------------|
| 014-20063 | Acatic Acid | 1 mL×5A |
| 018-20061 | Acetic Acid | 50 mL |
| 063-04533 | | 1 mL×5A |
| 067-04531 | Formic Acid (abt. 99%) | 50 mL |

for high-resolution MS (for QTofMS)

Absorbance

Gradient test*

Fluorescence test

Particles

Suitability for QTofMS analysis

*The gradient test for QTofMS was conducted as a QTofMS-compatible test.

Since a high-performance liquid chromatograph-mass spectrometer (LC-MS) has excellent characteristics, including the width of the scope of measuring objects and selectivity, it has become popular in a wide range of fields. Mass spectrometers with higher sensitivity and resolution have become popularized these days, and they are applied in research of metabolomics, proteomics, etc.




This product is a high-purity solvent optimal for measurements using high-sensitivity and high-resolution MS

✓ Testing requirement for compatibility with QTofMS

- A wide mass range by full scan is guaranteed. (mass range: 50–3000)
- Quality assurance with utilization of multivariate analysis-principal component analysis.
- ✓ Gradient measurement by column connection (UHPLC measurement)
- ✓ Particle assurance
- ✓ Content and container that ensure that the solvent can be used in high-purity conditions at any time
 - 500 mL or 1 L content optimal for using up
 - Glass bottles with special treatment that suppresses the elution of impurities are used (Patent No. WO/2005/047884)

Low background noise

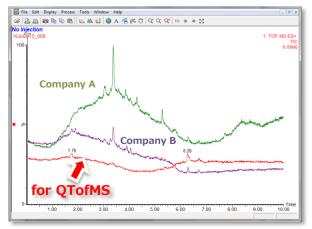

Comparison of background noise between our solvent for QTofMS and LC/MS solvents of other companies (Companies A and B)

We reviewed the manufacturing processes from various points so that it can be used in high-performance MS.

This product showed a lower background noise compared to the LC/MS solvents of other companies (Companies A and B) in quality confirmation using QTofMS.

[Total ion chromatogram]

Green : Company A Purple : Company B

Red : for QTofMS (Fujifilm Wako)

Small packaging for usability

Once opened, solvents become contaminated due to entry of air in laboratory, etc. To prevent contamination after opening, this product comes in 500 mL or 1 L content that can be used up.



These are the chromatograms immediately after opening (green) and 1 week after opening (red). Since it is highly likely that the increase in background noise would cause adverse effects on the measurement results in high-performance mass spectrometers, it is recommended that the mobile phase is used up each time.

Product list (for QTofMS)

| Product Number | Product Name | Package Size |
|----------------|-----------------|--------------|
| 018-26225 | Acatamituila | 500 mL |
| 016-26221 | Acetonitrile | 1 L |
| 130-18545 | Mothernel | 500 mL |
| 138-18541 | Methanol | 1 L |
| 164-27515 | 2-Propanol | 500 mL |
| 212-01601 | Ultrapure Water | 1 L |

Conditions and equipment for data acquisition

| Conditions for data acquisition | | |
|---------------------------------|--|--|
| HPLC | ACQUITY UPLC H-Class (Waters) | |
| Column | Wakopak [®] Ultra C18-2, 2.1mmφ×50mm | |
| Column temperature | 40℃±0.3℃ | |
| Eluate | Acetonitrile : Water [Gradient] | |
| Detector | UV 210nm | |
| Mass Spectrometer | Time-of-Flight Mass Spectrometer (Waters Xevo G2-XS) | |
| Ionization | ESI | |
| Mass Range | <i>m</i> / <i>z</i> 50∼3000 | |
| Polarity | positive | |



Waters Xevo G2-XS

for GPC

Absorbance

Gradient test

Suitability for GPC analysis

GPC solvents are solvents for which it has been confirmed that there is no coloring under heating for polymer dissolution or analysis. They are also guaranteed moisture, peroxides, non-volatile substances or the less changes in refractive index or ultraviolet absorption caused by impurities. This grades are optimal for preparation of GPC eluents.

Example of product specifications (1,1,1,3,3,3-Hexafluoro-2-propanol)

| REQUIREMENT | | SPECIFICATION |
|---|------------|------------------------|
| Appearance | | Colorless clear liquid |
| Density (20℃ | (g/mL) | 1.617-1.624 |
| Absorbance | 200 nm | max. 0.05 |
| | 220 nm | max. 0.03 |
| | 240 nm | max. 0.02 |
| | 254-400 nm | max. 0.01 |
| Suitability for | GPC test | to pass test |
| Water (%) | | max. 0.03 |
| Residue after evaporation (%) | | max. 0.001 |
| Acidity (as HF) (%) | | max. 0.001 |
| Peroxides (as H ₂ O ₂) (ppm) | | max. 5 |
| Assay (cGC) (%) | | min. 99.5 |

Product list (for GPC)

| Product Number | Product Name | Package Size |
|----------------|-----------------------------------|--------------|
| 036-24481 | Chlavafavor | 1 L |
| 032-24483 | Chloroform | 3 L |
| 039-24493 | Chloroform, Amylene added | 3 L |
| 043-33841 | o-Dichlorobenzene | 1 L |
| 049-33843 | 0-Dichiorobenzene | 3 L |
| 045-33921 | A/ A/ Dissable de ceta maido | 1 L |
| 041-33923 | N,N-Dimethylacetamide | 3 L |
| 046-33831 | N/ N/ Dissorby uforms provide | 1 L |
| 042-33833 | N,N-Dimethylformamide | 3 L |
| 048-33911 | Dimethyl Sulfoxide | 1 L |
| 044-33913 | | 3 L |
| 082-10311 | 1 1 1 2 2 2 110 (25) (25) | 100 mL |
| 084-10315 | 1,1,1,3,3,3-Hexafluoro-2-propanol | 500 mL |
| 134-18521 | 1-Methyl-2-pyrrolidone | 1 L |
| 130-18523 | | 3 L |
| 205-20071 | Totunhy duefy was Ctabilizer Free | 1 L |
| 201-20073 | Tetrahydrofuran, Stabilizer Free | 3 L |
| 209-20091 | T | 1 L |
| 205-20093 | Tetrahydrofuran, with Stabilizer | 3 L |

for the Japanese Pharmacopoeia General Tests


JP method

USP-NF method

EP method

The solvents for the Japanese Pharmacopoeia general tests (liquid chromatography) are HPLC solvents for which compliance with the Japanese Pharmacopoeia (JP) reagent standards, the United States Pharmacopoeia (USP) reagent standards, and the European Pharmacopoeia (EP) reagent standards is guaranteed, and for which some testing requirement of our company's own have been added.

Example of product label

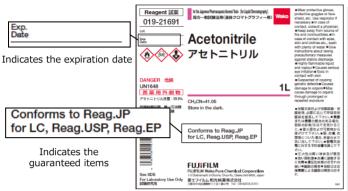

Example of product specifications (Acetonitrile)

| REQUIREMENT | JP method (*added requirement by us) | USP-NF method | EP method |
|--|--|------------------------|--------------------------------|
| Description | Colorless clear liquid | Colorless clear liquid | Colorless clear liquid |
| Color (APHA) | - | max. 10 | - |
| Density | 0.780 - 0.783 g/ml* | - | Sp. Gr. (20/20℃) about 0.78 |
| Refractive n_D^{20} | 1.343~1.346* | - | about 1.344 |
| Absorbance | - | - | - |
| 190 nm | - | max. 1.00 | - |
| 200 nm | max. 0.07 (0.05*) | - | max. 0.10 |
| 210 nm | max. 0.046 (0.03*) | - | - |
| 220 nm | max. 0.027 (0.02*) | max. 0.05 | - |
| 230 nm | max. 0.014 (0.01*) | - | - |
| 240 nm | max. 0.009 (0.005*) | - | - |
| 254 nm | - | max. 0.01 | - |
| 250-280 nm | - | max. 0.01 | - |
| 255-420 nm | - | - | max. 0.01 |
| 240-800 nm | - | - | max. 0.01 |
| Gradient test | to pass test* | to pass test | - |
| Residue after evaporation | 0.001 %* | 0.005%以下 | - |
| Titrable acid | Acidity 0.001 %* | 8 μeq/g以下 | - |
| Titrable base | - | 0.6 μeg/g以下 | - |
| Water | max. 0.05 %* | 0.003 | - |
| pH (100g/L) | - | - | Neutral to litmus paper |
| Distillation range (80-82 ℃) | - | - | min. 95 % |
| Ammonium (NH ₄) | 0.3 ppm* | - | - |
| Peroxides (H ₂ O ₂) | 5 ppm* | - | - |
| Substances reducing permanganate | to pass test* | - | - |
| Fluorescence test | to pass test* | - | - |
| Assay | min. 99.8 %* | min. 99.5 % | min. 99.9 % |

Product list (Solvents for the Japanese Pharmacopoeia General Tests)

| Product Number | Product Name | Package Size |
|----------------|----------------------------------|--------------|
| 019-21691 | | 1 L |
| 015-21693 | Acetonitrile | 3 L |
| 019-21696 | | 3 L×2 |
| 085-08711 | Hexane | 1 L |
| 081-08713 | | 3 L |
| 136-15661 | | 1 L |
| 132-15663 | Methanol | 3 L |
| 136-15666 | | 3 L×2 |
| 206-20981 | Tetrahydrofuran, Stabilizer Free | 1 L |

Protective Jacket Dedicated to Gallon Glass bottle GalloTect™

GalloTect™ is a protective jacket dedicated to gallon glass bottle. FUJIFILM Wako has developed GalloTect™ with the desire to add safety to daily testing and research activities, based on user feedback such as "when lifted the gallon glass bottle, it smashed and cracked."



- √ Carry out solvent resistance test with 11 kinds of solvents!
 - Solvent resistance of body material has been confirmed for a week with various solvents.
 - Tested solvents: Acetone, Acetonitrile, Benzyl Alcohol, 1-Butanol, DMF, DMSO, Ethanol, Glycerin, Methanol, 2-Propanol, Water
- ✓ A durability test was performed with a 10 kg weight! Weight test was performed at 10 kg for 1 hour.
- ✓ The lid was opened and closed 200 times to confirm that there is no loosening!
- ✓ All-plastics not requiring waste separation! A resin rivet was adopted for the screw portion. When discarding, a troublesome separation is not necessary!
- √ Just fit your reagent shelf! We focused on the design and eliminated the bulkiness thoroughly. The conventional reagent shelf can be stored gallon glass bottle with GalloTect™.













How to use

1. Place the gallon bottle and set the lid.



2. Turn the lid in the direction of the arrow until you hear a click.



3. Compleate



| Product Number | Product Name | Package Size |
|----------------|--------------|--------------|
| 293-36321 | Gallo Tect™ | 1 |

Listed products are intended for laboratory research use only, and not to be used for drug, food or human use. / Please visit FUJIFILM Wako Laboratory Chemicals site: https://labchem-wako.fujifilm.com/ / This leaflet may contain products that cannot be exported to your country due to regulations. / Bulk quote requests for some products are welcomed. Please contact us.

FUJIFILM Wako Laboratory Chemicals site https://labchem-wako.fujifilm.com

1-2, Doshomachi 3-Chome, Chuo-ku, Osaka 540-8605, Japan Tel: +81 6 6203 3741 Fax: +81 6 6203 1999 ffwk-cservise@fujifilm.com

FUJIFILM Wako Chemicals U.S.A. Corporation 1600 Bellwood Road, Richmond, VA 23237, U.S.A Toll-Free (U.S. only): +1 877 714 1920 Tel: +1 804 271 7677 Fax: +1 804 271 7791

wkuslabchem@fujifilm.com

FUJIFILM Wako Chemicals (Hong Kong) Limited FUJIFILM Wako (Guangzhou) Trading Corporation Room 1111, 11/F, International Trade Centre, 11-19 Sha Tsui Road,

Tsuen Wan, N.T., Hong Kong Tel: +852-2799-9019 Fax: +852-2799-9808 Tel: +852-2799-9019 F wkhk.info@fujifilm.com

FUJIFILM Wako Chemicals Europe GmbH

Fuggerstr 12, 41468 Neuss, Germany Tel: +49 2131 311 0 Fax: +49 2131 311 100 labchem_wkeu@fujifilm.com

Room 3003, 30/F., Dong Shan Plaza 69, Xian Lie Zhong Road Guangzhou, 510095, China Tel: +86-20-8732-6381(Guangzhou) Tel: +86-21-6288-4751(Shanghai) Tel: +86-10-6413-6388(Beijing) wkgz.info@fujifilm.com