







Anti Iba1, Goat Polyclonal Antibody

Application data -immunohistochemistry-

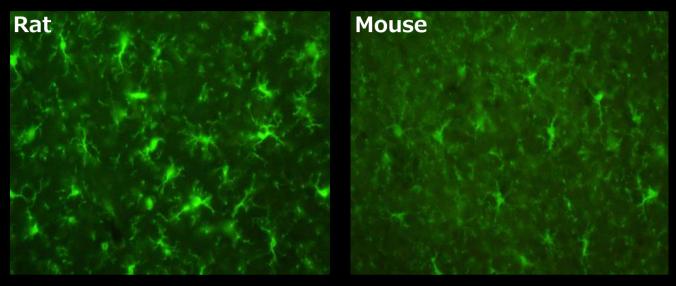

Immunohistochemistry of frozen section in rat(left) and mouse(right) brain cerebral cortex using anti Iba1, goat polyclonal antibody at dilution of 1:250, and anti goat IgG, Alexa Flour 488-conjugated. The data were provided by Dr. Nakajima from Soka University in Japan.

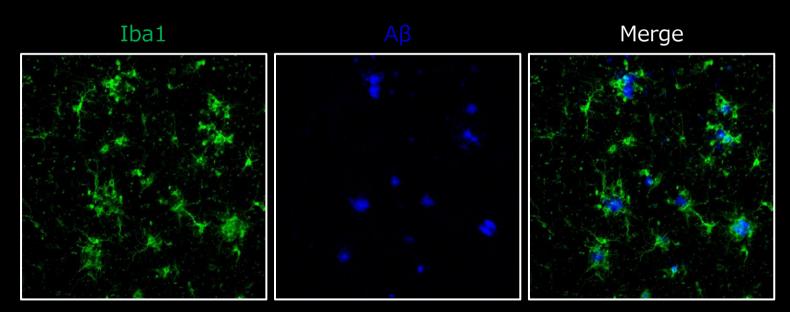

Immunohistochemistry of frozen section in mouse model of Alzheimer's disease (APP^{NL-G-F} mouse) brain cerebral neocortex using anti Iba1, goat polyclonal antibody at dilution of 1:1,000, and anti goat IgG, Alexa Flour 488-conjugated. A β were stained by 0.001% FSB solution(Amyloid fluorescent probe). The data were provided by Dr. Sakakibara, National Center for Geriatrics and Gerontology in Japan.

Product Name	Wako Cat. No	Pkg. Size	Grade	Storage Condition
Anti Iba1, Goat	011-27991	100µL	For Immunochemistry	Keep at -20°C

Application data -Western Blotting-

	Microglia	Neuron	Astrocyte	Cortex
→				

■ Sample

Rat primary microglia 10µgRat primary Neuron 10µg

· Rat primary Astrocyte 10μg

• Rat Brain Cortex 100µg

■ Antibody concentration 1/1000

The data were provided by Dr. Nakajima from Soka University in Japan.

Antibody Profile					
Antigen	Synthetic peptide corresponding to C-terminal of Iba1				
Buffer	TBS				
Host	Goat				
Species cross-reactivity	Rat and Mouse				
Antibody concentration	0.5-0.7 mg/mL				
Applications	Immunohistochemistry(frozen section) 1:250-1,000 Immunohistochemistry(paraffin section) 1:250-1,000 Western blotting 1:1,000				

Product List

Product Name	Wako Cat. No	Pkg. Size	Grade	Storage Condition
Anti Iba1, Rabbit (for Immunocytochemistry)	019-19741	50µg	For Immunochemistry	Keep at -20℃
Anti Iba1, Rabbit(for Paraffin Section)	013-27691	50µg	For Immunochemistry	Keep at -20℃
Anti Iba1, Rabbit, Biotin-conjugated	016-26461	100µL	For Immunochemistry	Keep at 2-10℃
Anti Iba1, Rabbit, Red Fluorochrome(635)-conjugated	013-26471	100µL	For Immunochemistry	Keep at 2-10℃
Anti Iba1, Rabbit (for Western Blotting)	016-20001	50µg	For Immunochemistry	Keep at -20℃
Anti Iba1, Monoclonal Antibody(NCNP24)	012-26723 016-26721	10μL 50μL	For Immunochemistry	Keep at -20℃
Anti Human Iba1, Monoclonal Antibody(NCNP27)	017-27591 013-27593	10μL 50μL	For Immunochemistry	Keep at -20℃

Listed products are intended for laboratory research use only, and not to be used for drug, food or human use. / Please visit our online catalog to search for other products from Wako; https://labchem-wako.fujifilm.com / This leaflet may contain products that cannot be exported to your country due to regulations. / Bulk quote requests for some products are welcomed. Please contact us.

FUJIFILM Wako Pure Chemical Corporation

https://labchem-wako.fujifilm.com 1-2, Doshomachi 3-Chome Chuo-Ku, Osaka 540-8605, Japan Tel.: +81 6 6203 3741 Fax: +81 6 6203 1999 ffwk-cservice@fujifilm.com

FUJIFILM Wako Chemicals U.S.A. Corporation

www.wakousa.com 1600 Bellwood Road Richmond, VA 23237, USA Toll-Free (U.S. only): $+1\,877\,714\,1920$ Tel.: $+1\,804\,714\,1920$ Fax: $+1\,804\,271\,7791$ wkuslabchem@fujiflim.com



FUJIFILM Wako Chemicals Europe GmbH

www.wako-chemicals.de Fuggerstr 12, 41468 Neuss, Germany Tel.: +49 2131 311 0 Fax: +49 2131 311 100 labchem_wkeu@fujifilm.com