

Description

- Origin :** Broadbean Peptone is obtained by enzymatic hydrolysis of broadbean protein concentrates.
- Context :** this peptone is animal free and is unconcerned by GMO issues (European Directive 2001/18/CE). It is, to the best of our knowledge, non-allergenic and not subject to mandatory declaration acc. to Annex IIIa of the EU directive 2003/89/EC, updated with 2006/142/EC. Product is **Kosher DE** and **Halal** certified. An identical **Kosher Pareve Halal** reference also available under reference **A223100**.
- Application :** This peptone demonstrates excellent growth of a wide variety of microorganisms, particularly lactic acid bacteria (LAB) and may serve as a substitute for meat-based peptones in certain circumstances, in both industrial fermentation media, diagnostic applications and in cell culture.



Physical properties

- Appearance :** beige powder
Stability (2% in solution) : stable
Solubility in water at 2% : total

Microbiological controls

Total aerobic mesophilic flora ≤ 5000 cfu/g

Chemical analysis

- Total nitrogen (N_T) :** 10.0%
α-amino nitrogen (N_α) : 2.8%
N_α / N_T : 0.28
Sulfuric ash : 12.6 %
pH (2% in solution) : 6.9
Total carbohydrates : approx. 11.7%
Chlorides (as NaCl) : 0.3%
Loss on drying : ≤ 6.0%


Standard packaging

25 kg carton ; other formats inquire.
 Delivered with Certificate of Analysis, Certificate of Origin. Separate non-Animal, GMO, Kosher, Halal and allergen statements available upon request.

Amino acid (AA) distribution (mg/g)

	Total AA	Free AA		Total AA	Free AA
Aspartic acid	73.0	2.24	Methionine	10.5	0.51
Threonine	24.0	2.71	Isoleucine	26.5	1.82
Serine	32.2	3.13	Leucine	26.5	7.62
Glutamic acid	117.0	6.09	Tyrosine	26.0	2.50
Proline	28.0	1.06	Phenylalanine	31.5	5.70
Glycine	26.0	0.64	Histidine	17.0	2.32
Alanine	30.5	2.17	Lysine	41.0	3.17
Cysteine	5.6	1.20	Arginine	63.5	18.06
Valine	29.0	4.57	Tryptophan	4.5	0.00

Molecular weight distribution (Daltons)

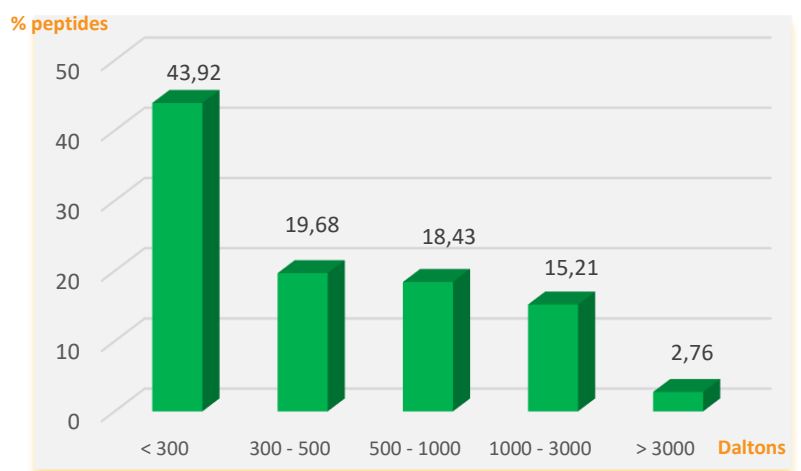



Sanitary / Regulatory statement

This plant peptone is classified animal-free by SOLABIA S.A.S. Based on the manufacturing protocol, we attest that no animal raw materials are prescribed for use in the production this product nor are any of the raw materials derived from animal products. Also, to the best of our knowledge, the product contains no genetically modified organisms as defined by current legislation for labelling (absence = less than 0.9%).

Storage


Keep in original packaging when not in use, in a dry area ideally between 10 and 35°C. Avoid direct sunlight. Hygroscopic product. Expiry date : 5 years from date of manufacture.

OBSERVED MICROBIAL GROWTH POTENTIAL :

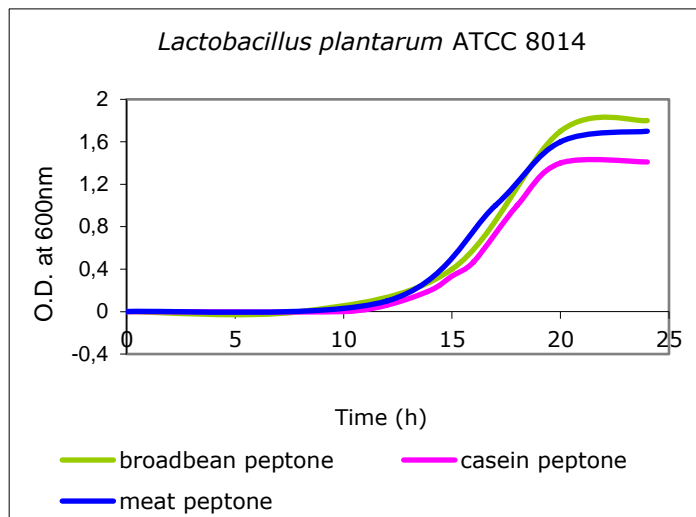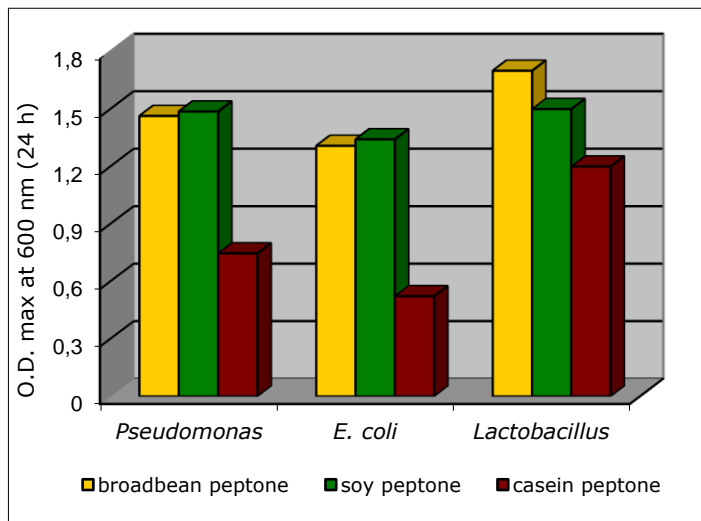

Inoculum 10^8 cfu / mL
 Growth medium: 3 % peptone + 0.25 % glucose
 (*Pseudomonas* & *E. coli*)
 Modified MRS medium, pH 6.4 (*Lactobacillus*)

Inoculum 10^2 cfu / mL
 Modified MRS medium, pH 6.4

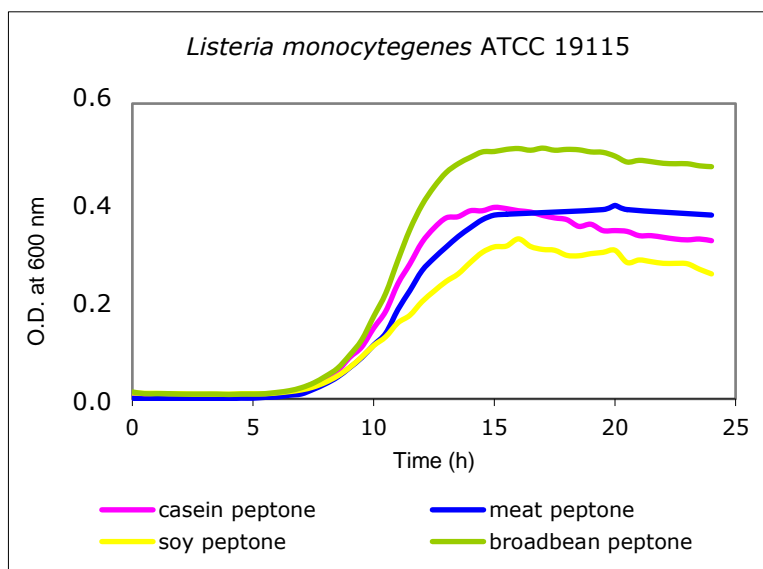

Inoculum 10^2 cfu/mL
 Culture medium : 1% peptone + 0.5% NaCl

Laboratory tests demonstrate excellent growth of lactic acid bacteria and *Enterobacteria*, used as either a replacement for bovine substrates or as a stand-alone peptone. Results may differ depending on individual laboratory conditions and for other genera, species and strains.



Certifications profile :

Produced under **ISO 9001 v 2015** certification



Manufacturing site and quality system open to audits by qualified customers*.

CoA available online : use product code **A220100** + lot number
<https://www.solabia.com/coa>

Allergen considerations : **NO** GMO issues : **NO**

Kosher certification : yes

Orthodox Union ; New York, NY USA.
 Rabbi Menachem Adler, Administrator

Halal certification : Yes A220100
Yes A223100



HFFIA ; The Hague, Netherlands.
 Inquire with Solabia for more information

V. 02/2023

The information presented in this document is submitted in good faith based on internal testing performed at Solabia S.A.S. and represents the best of our knowledge at the present time. It is provided as a guide and no warranty, implied or otherwise is associated with this data, nor is any liability assumed for patent infringement. All data represents typical analyses not to be taken for exact specifications. *Audits may involve fees dependent on volume of business.

End-users are directed to perform proprietary tests to determine suitability and performance for specific applications. The information and results contained in this technical data sheet are susceptible to modification at any time, without warning.