

Special Section – New Models in Drug Metabolism and Transport

A Novel In Vitro Experimental System for the Evaluation of Drug Metabolism: Cofactor-Supplemented Permeabilized Cryopreserved Human Hepatocytes (MetMax Cryopreserved Human Hepatocytes)[§]

Albert P. Li, Ming-Chih David Ho, Kirsten Amaral, and Carol Loretz

In Vitro ADMET Laboratories Inc., Columbia, Maryland (A.P.L., C.L.) and In Vitro ADMET Laboratories Inc., Malden, Massachusetts (M.-C.D.H., K.A.)

Received November 22, 2017; accepted January 22, 2018

ABSTRACT

We report here a novel experimental system, cryopreserved MetMax human hepatocytes (MMHHs), for in vitro drug metabolism studies. MMHHs consist of cofactor-supplemented permeabilized cryopreserved human hepatocytes. The use procedures for MMHHs are significantly simplified from that for conventional cryopreserved human hepatocytes (CCHHs): 1) storage at -80°C instead of in liquid nitrogen and 2) usage directly after thawing without centrifugation and microscopic evaluation of cell density and viability and cell density adjustment. In this study, we compared MMHHs and CCHHs in CYP1A2, CYP2A6, CYP2B6, CYP2C8, CYP2C9, CYP2D6, CYP2E1, CYP3A4, CYP2J2, monoamine oxidase A, aldehyde oxidase, flavin-containing monooxygenase, UDP-glucuronyl transferase, *SULT*, *N*-acetyltransferase 1, and acetaminophen glutathione (GSH) conjugation activities based on liquid chromatography–tandem mass spectrometry quantification of substrate metabolism. MMHHs were prepared from CCHHs consisting of hepatocytes pooled from

10 individual donors. The drug metabolizing enzyme activities of both CCHHs and MMHHs were cell concentration and time dependent, with specific activities of MMHHs ranging from 27.2% (carboxylesterase 2) to 234.2% (acetaminophen GSH conjugation) of that for CCHHs. As observed in CCHHs, sequential oxidation and conjugation was observed in MMHHs for coumarin, 7-ethoxycoumarin, and acetaminophen. 7-Hydroxycoumarin conjugation results showed that metabolic pathways in MMHHs could be selected via the choice of cofactors, with glucuronidation but not sulfation observed in the presence of UDP-glucuronic acid and not 3-phosphoadenosine-5-phosphosulfate, and vice versa. Results with noncytotoxic and cytotoxic concentrations of acetaminophen showed that drug metabolism was compromised in CCHHs but not in MMHHs. Our results suggest that the MMHHs system represents a convenient and robust in vitro experimental system for the evaluation of drug metabolism.

Introduction

Successful cryopreservation of human hepatocytes is a major reason for the routine application of this experimental system for preclinical assessment of human-specific drug properties including metabolic fate, drug-drug interactions, and drug toxicity (Li, 2007; Godoy et al., 2013; Jones et al., 2017b; Wood et al., 2017). Since the functional performance of cryopreserved human hepatocytes can be significantly compromised by cellular damage, use of cryopreserved human hepatocytes requires practitioners to have extensive experience with handling procedures, including thawing, centrifugation, resuspension, microscopic quantification of cell viability, adjustment of cell concentration, and the ultimate delivery of the cells into experimental vessels for experimentation. Furthermore, the use of cryopreserved hepatocytes requires liquid nitrogen cryogenic freezers for storage, which may

not be readily available in laboratories where in vitro drug metabolism studies are routinely performed.

Before the development of cryopreserved hepatocyte technologies, hepatic subcellular fractions such as human liver microsomes and postmitochondrial supernatants (S9 or S10), were used exclusively in drug metabolism studies (Iwatsubo et al., 1997; Gombar et al., 2003). These subcellular systems, especially human liver microsomes, continue to be practical and useful experimental systems in pharmaceutical industrial laboratories, including screening of new chemical entities for metabolic stability (Halladay et al., 2007; Choi et al., 2015), estimation of in vivo hepatic clearance (Obach, 2011; Chen et al., 2017), evaluation of cytochrome P450–related drug properties (Dinger et al., 2014), and UDPGT–mediated drug metabolism (Walsky et al., 2012; Joo et al., 2014). Compared with cryopreserved hepatocytes, the use of subcellular fractions requires relatively simple application procedures, and they are relatively robust and are not readily subjected to functional damages due to handling. A major drawback of the use of subcellular fractions is the incompleteness

<https://doi.org/10.1124/dmd.117.079657>

[§]This article has supplemental material available at dmd.aspetjournals.org.

ABBREVIATIONS: CCHH, conventional cryopreserved human hepatocyte; DME, drug metabolizing enzyme; LC/MS-MS, liquid chromatography–tandem mass spectrometry; MMHH, MetMax human hepatocyte; PAPS, 3-phosphoadenosine-5-phosphosulfate; *SULT*, sulfotransferase; UDPGA, UDP-glucuronic acid.

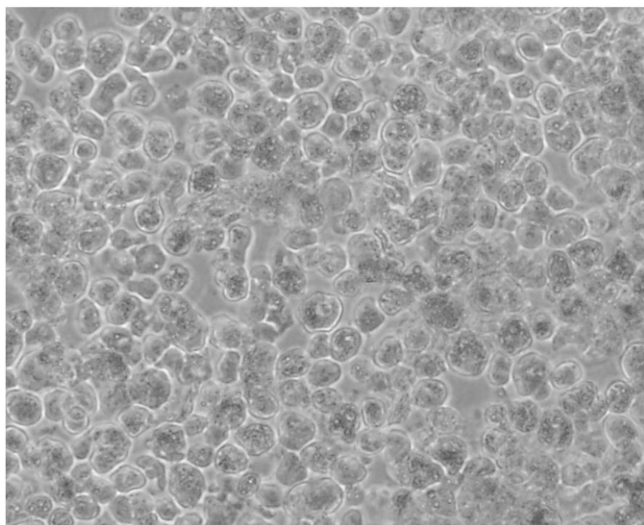


Fig. 1. Phase-contrast microscopy of MMHs. The hepatocytes were plated into a 24-well plate and allowed to settle for approximately 10 minutes. The cells exhibited intact plasma membranes similar to CCHHs.

of the drug metabolizing enzyme (DME) pathways, thus the data obtained may not allow accurate assessment of *in vivo* hepatic metabolic fate. In contrast, cryopreserved hepatocytes contain complete, undisrupted DMEs and cofactors, and are considered to be the gold standard for the evaluation of *in vitro* human drug metabolism (Fabre et al., 1990; Ulrich et al., 1995; Hewitt et al., 2007). For this reason, despite the relatively more complicated experimental procedures, cryopreserved human hepatocytes have been applied in studies routinely performed with subcellular fractions, including hepatic clearance (Di et al., 2013; Peng et al., 2016) and cytochrome P450 inhibition (Doshi and Li, 2011; Li and Doshi, 2011; Kazmi et al., 2015), to obtain data that are more likely to reflect *in vivo* events.

It would be desirable to have an experimental system that combines the completeness of the DME pathways in hepatocytes, and the robustness and simplicity of the application procedures of subcellular fractions. Toward this goal we have developed a novel experimental system that uses cryopreserved MetMax human hepatocytes (MMHs). MMHs are permeabilized cryopreserved human hepatocytes supplemented with key DME cofactors. Similar to subcellular fractions such as human liver microsomes, MMHs can be stored at -80°C in contrast to conventional cryopreserved human hepatocytes (CCHHs) which require liquid nitrogen storage. Furthermore, similar to subcellular fractions, MMHs can be used directly after thawing without centrifugation and microscopic evaluation of viability and cell concentration as required for CCHHs.

We report here a comparison of MMHs and CCHHs in phase 1 oxidation and phase 2 conjugation DME activities as part of our ongoing research program to evaluate the applicability of this novel experimental system to study drug metabolism.

Materials and Methods

Chemicals. Dextrophan tartrate, diclofenac sodium salt, 4-hydroxydiclofenac, S-mephenytoin, 4-hydroxyquinoline, paclitaxel, and testosterone were purchased from Cayman Chemical (Ann Arbor, MI). 7-Hydroxycoumarin was purchased from Chem Service (West Chester, PA). Benzydamine *N*-oxide, 7-hydroxycoumarin sulfate potassium salt, kynuramine hydrobromide, and *N*-acetyl sulfamethazine were obtained from Santa Cruz Biotechnology (Dallas, TX). 4-Acetamidobenzoic acid, *p*-acetamidophenyl β -D-glucuronide sodium salt, 4-aminobenzoic acid, benzydamine hydrochloride, chlorzoxazone, coumarin, dextromethorphan hydrobromide, 6 β -hydroxytestosterone, 7-hydroxycoumarin β -D-glucuronide sodium salt, 7-ethoxycoumarin, paracetamol sulfate potassium, phenacetin, and sulfamethazine were purchased from Sigma Aldrich (St. Louis, MO). Carbazeran, 4-hydroxycarbazeran, 6-hydroxychlorzoxazone, 6 α -hydroxypaclitaxel, acetaminophen glutathione disodium salt, midazolam, 1'-hydroxymidazolam, and 4-hydroxy-S-mephenytoin were obtained from Toronto Research Chemicals (Toronto, Canada). All other DME substrates were obtained from Sigma Aldrich.

CCHH and MMHH Preparations. Cryopreserved human hepatocytes previously prepared from 10 individual donors were used in the preparation of

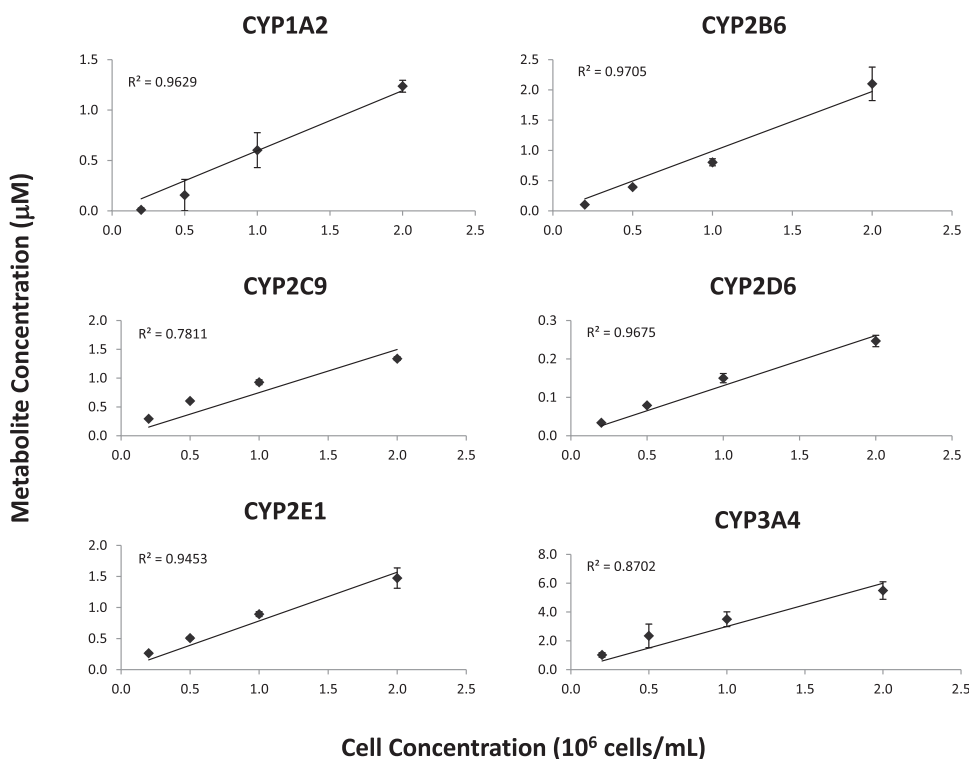


Fig. 2. Metabolite formation as a function of cell concentration in MMHs. Cytochrome P450 isoform-selective pathways evaluated were CYP1A2 (phenacetin *O*-deethylation), CYP2B6 (bupropion hydroxylation), CYP2C9 (diclofenac 4'-hydroxylation), CYP2D6 (dextromethorphan *O*-demethylation), CYP2E1 (chlorzoxazone 6-hydroxylation), and CYP3A4 (midazolam 1'-hydroxylation). The cell concentrations evaluated were 0.25, 0.5, 1, and 2×10^6 cells/mL. (Error bars: S.D. of triplicate determinations).

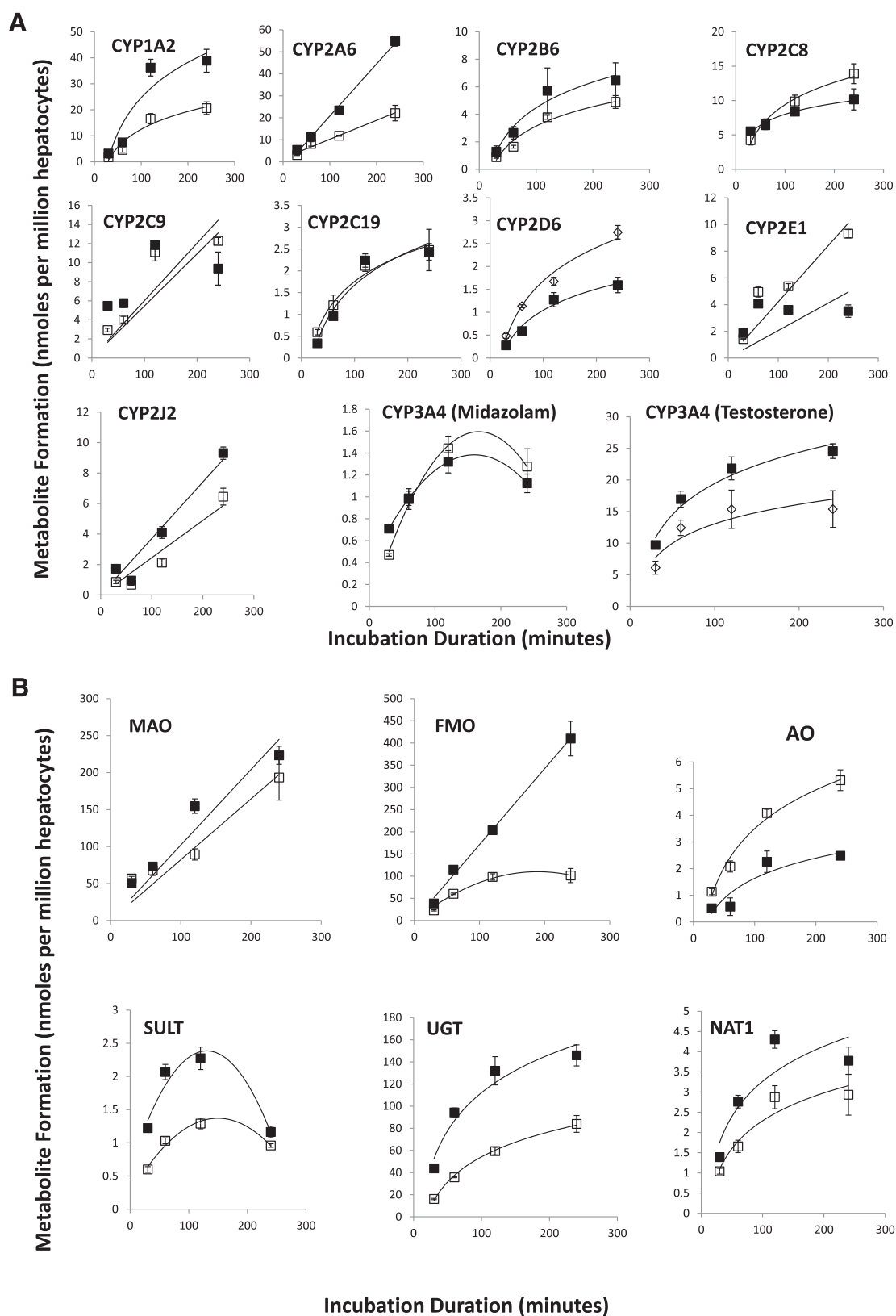


Fig. 3. Metabolite formation as a function of incubation durations in CCHHs (open squares) and MMHs (filled squares). The time durations evaluated were 30, 60, 120, and 240 minutes. The cell concentration used was 1×10^6 cells/ml. Cytochrome P450 (A) and non-cytochrome P450 (B) pathway-selective substrate metabolism are shown. (Error bars: S.D. of triplicate determinations).

the pooled donor CCHHs and MMHs. Pooled CCHHs were prepared by thawing of cryopreserved human hepatocytes from individual lots, pooling of the individual lots, and recryopreserving the pooled multiple donor hepatocytes at a

cell density of approximately 5×10^6 cells/ml using a patented technology (Li, 2015a). For the preparation of CCHHs, hepatocytes from each donor were isolated via collagen digestion of human livers and cryopreserved in liquid

Conventional Cryopreserved Human Hepatocytes


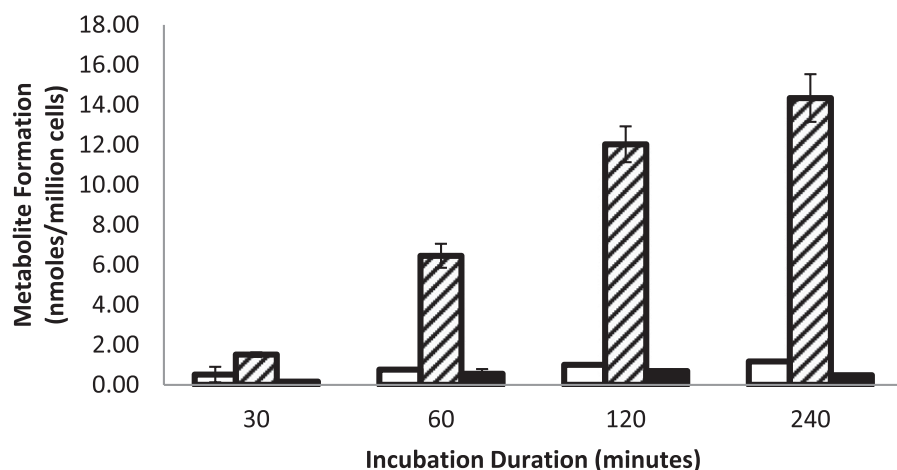


Fig. 4. A comparison of CCHHs (top) and MMHHs (bottom) in the formation of secondary metabolites from 7-ethoxycoumarin. The metabolites quantified were 7-hydroxycoumarin (7HC; open bars) and its glucuronide (7HCG; shaded bars) and sulfate (7HCS; filled bars) conjugates. The incubation durations were 30, 60, 120, and 240 minutes. The cell concentration used was 1×10^6 cells/ml. (Error bars: S.D. of triplicate determinations).

MetMax Cryopreserved Human Hepatocytes

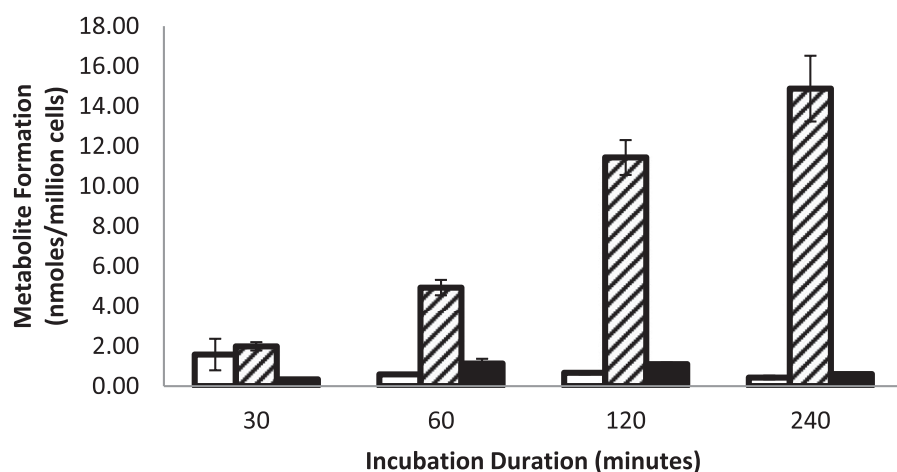

nitrogen (Li, 2017). Pooled MMHHs were prepared from the same lot of CCHHs used in the study via thawing and recovery of the CCHH lot, recovery of the hepatocytes followed by permeabilization of the cells with a proprietary technology (patent pending), and re-cryopreserved at a cell

density of 2×10^6 cells/ml. The human livers used for hepatocyte isolation were obtained from the International Institute for the Advancement of Medicine (Edison, NJ) and the National Disease Research Exchange (Philadelphia, PA). The demographics of the 10 donors are given in Supplemental Table 1.

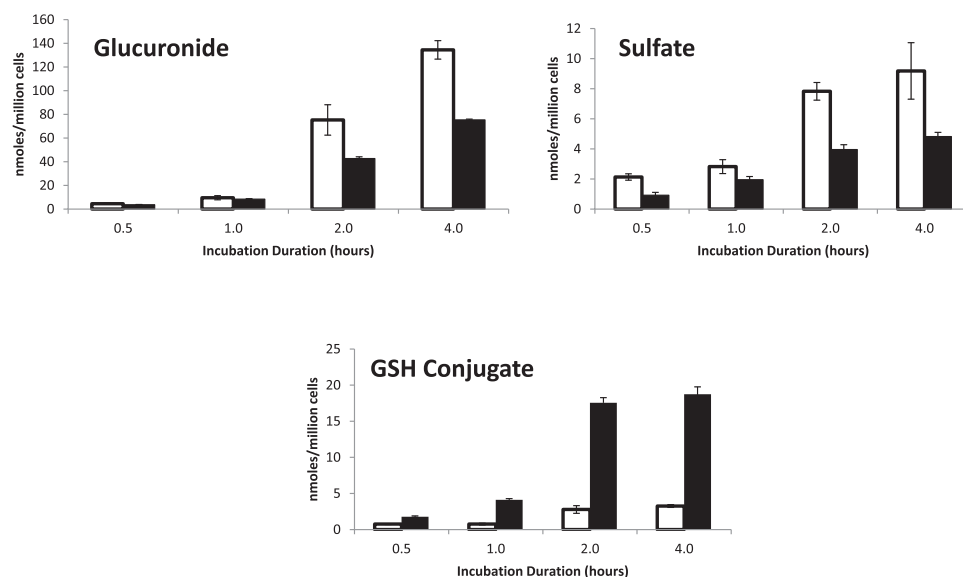


Fig. 5. Time-dependent formation of acetaminophen conjugates in CCHHs (open bars) and MMHHs (filled bars) at the noncytotoxic concentration of 10 mM. The metabolites quantified were glucuronide (A), sulfate (B), and glutathione (GSH) (C) conjugates. (Error bars: S.D. of triplicate determinations).

TABLE 1

Comparison of CCHHs in the metabolism of acetaminophen at the noncytotoxic concentration of 10 mM
Mean and S.D. of triplicate samples are shown.

| Incubation Time | Total Metabolite Formation in CCHHs (nmoles per Million Hepatocytes) | | | | | |
|-----------------|--|-------|--------------|------|------------------|------|
| | APAP Glucuronide | | APAP Sulfate | | APAP Glutathione | |
| | Mean | S.D. | Mean | S.D. | Mean | S.D. |
| <i>min</i> | | | | | | |
| 30 | 4.51 | 0.46 | 2.13 | 0.21 | 0.76 | 0.06 |
| 60 | 9.57 | 1.76 | 2.82 | 0.46 | 0.76 | 0.12 |
| 120 | 75.26 | 12.84 | 7.83 | 0.59 | 2.78 | 0.53 |
| 240 | 134.46 | 7.78 | 9.18 | 1.88 | 3.26 | 0.20 |

Incubation with DME Substrates. All incubations were performed in 96-well cell culture plates (Falcon; obtained from VWR Inc., Bridgeport, NJ). For CCHHs, the hepatocytes were thawed in a 37°C water bath, with 1 ml of the thawed suspension added to 50 ml of universal cryopreservation recovery medium (In Vitro ADMET Laboratories Inc., Columbia, MD) in a 50 ml conical cell culture tube (Falcon; obtained from VWR Inc.) and centrifuged at 100g for 10 minutes. The cell pellet from each conical tube was resuspended in 4 ml of HQM Hepatocyte Incubation Medium (In Vitro ADMET Laboratories Inc.) for viability determination (Trypan blue exclusion) and cell concentration determination. The cell suspension was then adjusted with HQM to 2× the final cell density. For MMHs, the hepatocytes were thawed and used directly without centrifugation or viability/cell concentration determination. All DME substrates were prepared in HQM at 2× the final concentrations and added at a volume of 50 µl per well in a 96-well cell culture plate. The CCHHs, MMHs, and substrate plates were prewarmed to 37°C for 15 minutes in a cell culture incubator before initiation of the incubation by pipetting 50 µl of CCHHs or MMHs into each well of the 96-well plates containing the substrates and returned to the cell culture incubator without shaking for the desired durations (30, 60, 90, 120, and 240 minutes). After incubation, 100 µl of acetonitrile was added into each well to terminate metabolism. After termination, the plates were stored at -80°C in a freezer for later liquid chromatography–tandem mass spectrometry (LC-MS/MS) quantification of metabolite formation. All DME activity studies were performed at a final cell density of 1×10^6 cells/ml except for the evaluation of cell density versus metabolite formation where CCHHs and MMHs were first prepared at cell concentrations of 4×10^6 cells per ml, followed by dilution to 2×10^6 , 1×10^6 , and 0.5×10^6 per ml before addition to the substrate plates, thereby yielding final concentrations of 2×10^6 , 1×10^6 , 0.5×10^6 , and 0.25×10^6 ml.

LC-MS/MS Analysis. Upon thawing, an aliquot of 200 µl of each sample was transferred from each well into a labeled 96-well plate followed by addition of 100 µl of acetonitrile solution containing the internal standard Tolbutamide (250 nM) and mixing. All samples were centrifuged at 3500 rpm for 5 minutes. An aliquot of 100 µl of supernatant from each was transferred to a 96-well plate and diluted with 200 µl of deionized water with mixing before LC-MS/MS analysis. CYP1A2, CYP2A6, CYP2B6, CYP2C8, CYP2C9, CYP2C19, CYP2D6, CYP2E1, CYP3A4 (midazolam 1'-hydroxylation), CYP3A4 (testosterone 6β-hydroxylation), 7-ethoxycoumarin O-deethylation, UDP-glucuronyl transferase, sulfotransferase (SULT), glutathione S-transferase, flavin-containing monooxygenase, monoamine oxidase, aldehyde oxidase and, N-acetyltransferase

1 metabolites were quantified by an API 4000 QTRAP mass spectrometer with an electrospray ionization source (AB SCIEX, Framingham, MA) connected to an Agilent 1200 Series HPLC (Agilent Technologies, Santa Clara, CA) using the LC-MS/MS multiple reactions monitoring mode and monitoring the mass transitions (parent-to-daughter ion) (Supplemental Table 2). A Zorbax Eclipse Plus C18 column (4.6 × 75 mm i.d., 3.5 µm; Agilent Technologies) at a flow rate of 1 ml/min was used for chromatography separation. The mobile phase consisted of 0.1% formic acid in acetonitrile (A) and 0.1% formic acid in water (B). The gradient for positive ion mode operation was programmed as: 0–2.5 minutes, increase B from 5% to 95%; 2.5–3.5 minutes, 95% B; 3.5–3.6 minutes, decrease B to 5%; run time, 5 minutes. The gradient program for the negative ion mode was: 0–3 minutes, increase B from 5% to 95%; 3–4 minutes, 95% B; 4–4.2 minutes, decrease B to 5%; run time 6 minutes. Data acquisition and data processing were performed with the software Analyst 1.6.2 (AB SCIEX).

Data Analysis. Data are presented as mean and S.D. of triplicate incubations derived using the Microsoft Excel 6.0 software. Statistical analysis was performed using the Student's *t* test with the Microsoft Excel 6.0 software, with the probability of *P* < 0.05 to be considered statistically significant. Specific activity (picomoles per minute per million hepatocytes) of each DME pathway was determined by dividing the total metabolite formed by the incubation time and normalized to cell concentration.

Results

Cell Morphology

When placed into a cell culture plate, the hepatocytes from the MMHs appeared as intact cells, which rapidly settled on the bottom of the wells. Phase contrast photomicroscopy showed that the MMHs were intact (not lysed) cells with distinct cell membranes (Fig. 1), similar to the CCHHs (data not shown). Trypan blue exclusion evaluation, however, showed that 100% of the cells would include the dye (data not shown), demonstrating that their plasma membranes were permeabilized.

Cell Concentration versus Metabolite Formation

MMHs were incubated for 30 minutes with various cytochrome P450 isoform-selective substrates at 0.25, 0.5, 1.0, and 2.0×10^6 cells/ml.

TABLE 2

Comparison of MMHs in the metabolism of acetaminophen at the noncytotoxic concentration of 10 mM
Mean and S.D. of triplicate samples are shown.

| Incubation Time | Total Metabolite Formation in MMHs (nmoles per Million Hepatocytes) | | | | | |
|-----------------|---|------|--------------|------|------------------|------|
| | APAP Glucuronide | | APAP Sulfate | | APAP Glutathione | |
| | Mean | S.D. | Mean | S.D. | Mean | S.D. |
| <i>min</i> | | | | | | |
| 30 | 3.62 | 0.29 | 0.92 | 0.19 | 1.78 | 0.10 |
| 60 | 8.40 | 0.41 | 1.95 | 0.21 | 4.12 | 0.17 |
| 120 | 42.66 | 1.46 | 3.95 | 0.32 | 17.56 | 0.70 |
| 240 | 75.32 | 0.68 | 4.82 | 0.28 | 18.74 | 1.01 |

TABLE 3

Comparison of CCHHs and MMHHs in acetaminophen metabolism at the lesser cytotoxic concentration of 100 mM and the cytotoxic concentration of 200 mM

Mean and S.D. of triplicate samples are shown.

| Substrate Concentration (μM) | Marker Metabolite | Total Metabolite Formation (pmol per Million Hepatocytes) | | | |
|------------------------------|---------------------------|---|-------|---------|-------|
| | | CCHH | | MMHH | |
| | | Average | S.D. | Average | S.D. |
| <i>mM</i> | | | | | |
| 100 | Acetaminophen Glutathione | 206.50 | 3.91 | 856.72 | 71.28 |
| 200 | Acetaminophen Glutathione | 99.81 | 7.49 | 676.70 | 54.43 |
| 100 | Acetaminophen Glucuronide | 572.55 | 27.12 | 1298.77 | 64.29 |
| 200 | Acetaminophen Glucuronide | 170.28 | 16.38 | 1361.84 | 83.36 |
| 100 | Acetaminophen Sulfate | 105.55 | 3.36 | 150.22 | 6.68 |
| 200 | Acetaminophen Sulfate | 52.58 | 1.40 | 156.36 | 11.86 |

Metabolite formation was found to increase with cell concentration for all of the substrates evaluated (Fig. 2). From the results, 1.0×10^6 cells/ml was chosen as the cell concentration for subsequent studies.

Incubation Duration versus Metabolite Formation

MMHHs and CCHHs were incubated with various pathway-selective substrates for 30, 60, 120, and 240 minutes. In general, similar relationships between metabolite formation and incubation duration were observed for the two hepatocyte systems, with linear or curve-linear time course (Fig. 3).

Secondary Metabolite Formation.

1. 7-Ethoxycoumarin metabolism: MMHHs and CCHHs were incubated with 7-ethoxycoumarin. 7-Ethoxycoumarin is first oxidized to 7-hydroxycoumarin, which in turn is subjected to metabolism by UDPGT and SULT to 7-hydroxycoumarin glucuronide and 7-hydroxycoumarin sulfate, respectively. Similar metabolite profiles were observed for MMHHs and CCHHs, with 7-hydroxycoumarin found to be extensively conjugated to the glucuronide and sulfate conjugates (Fig. 4).
2. Acetaminophen metabolism: Acetaminophen metabolism by MMHHs and CCHHs was evaluated at the noncytotoxic concentration of 10 mM at incubation durations of 30, 60, 120, and 240 minutes (Fig. 5; Tables 1 and 2). Time-dependent formation of glucuronide, sulfate, and glutathione conjugates was

observed in both MMHHs and CCHHs. The glutathione conjugate of acetaminophen formed in MMHHs ranged from approximately 200% at 30 minutes to approximately 600% of that formed in CCHHs at 240 minutes. In a separate study, MMHHs and CCHHs were incubated with 100 and 200 mM of acetaminophen (Table 3). The results confirmed the observation made with 10 mM acetaminophen, in that higher levels of acetaminophen glutathione were formed in MMHHs. Significant reduction of acetaminophen metabolism was observed at 200 mM in CCHHs, which was contributed to the cytotoxic effects at this drug concentration. No apparent reduction in metabolite formation was observed in MMHHs at 200 mM of acetaminophen.

Cofactor-Mediated Selection of Metabolic Pathways in MMHHs

MMHHs were prepared without cofactors and then incubated with coumarin with the following cofactor compositions: 1) NADPH [without UDP-glucuronic acid (UDPGA) or 3-phosphoadenosine-5-phosphosulfate (PAPS)]; 2) NADPH plus UDPGA (without PAPS); and 3) NADPH plus PAPS (without UDPGA). 7-Hydroxycoumarin was formed under all three cofactor compositions. 7-Hydroxy and 7-glucuronide (but not the 7-sulfate metabolites) were observed under condition 2; 7-hydroxylation and 7-sulfate (but not the 7-glucuronide metabolites) were observed under condition 3 (Fig. 6).

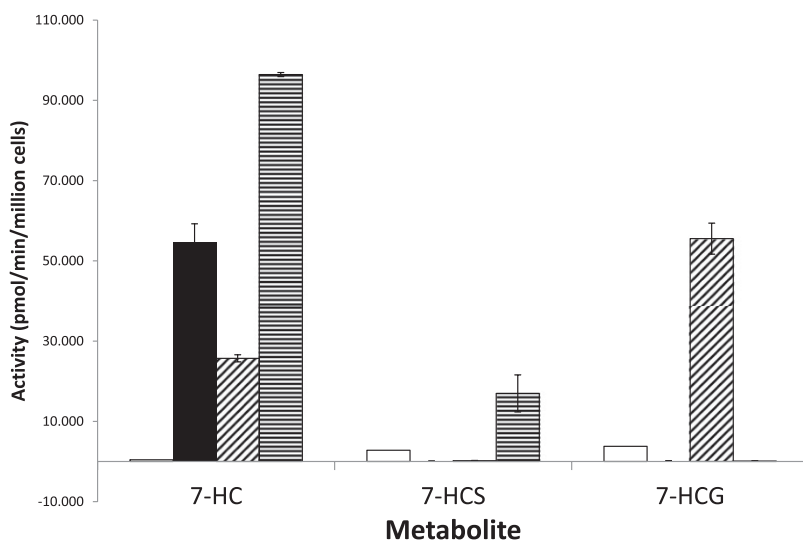


Fig. 6. Effects of cofactor composition on metabolite formation from coumarin in MMHHs. MMHHs were incubated with coumarin under four cofactor conditions: 1) no cofactor supplements (open bars), 2) NADPH only (filled bars), 3) NADPH and UDPGA (slanted shaded bars), and 4) NADPH and PAPS (horizontally shaded bars). (Error bars: S.D. of triplicate determinations).

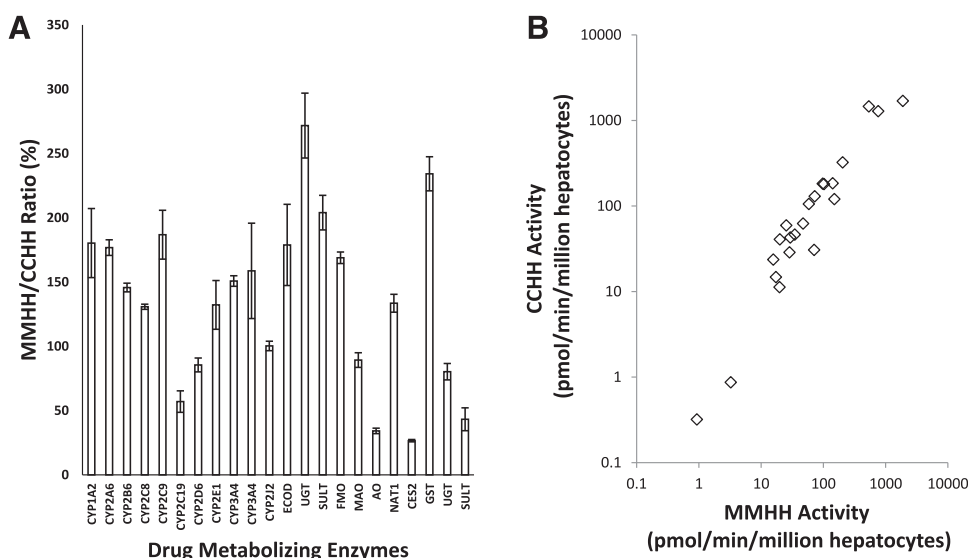


Fig. 7. A comparison of CCHHs and MMHHs in DME activities: DME activities of MMHHs expressed as percentages of that for CCHHs (A), calculated using the following equation: $[\text{Activity (MMHH)}/\text{Activity (CCHH)}] \times 100\%$, and correlation of MMHH and CCHH activities (B). Each data point in (B) represents activity data for each of the DMEs identified in (A). Numerical data are presented in Supplemental Table 3. (Error bars: S.D. of triplicate determinations).

Activity Comparison between MMHHs and CCHHs

A comparison of the rates of metabolism of the various drug metabolism enzyme-selective substrates in MMHHs is shown in Fig. 7, with MMHHs activities expressed as percentages of that for CCHHs (Fig. 7A) and MMHH activities plotted versus CCHH activities (Fig. 7B). Numerical data for DME activities for MMHHs and CCHHs are presented in Supplemental Table 3. The results demonstrate that MMHHs were active in all of the enzyme pathways evaluated, with activities in general similar to that observed in CCHHs.

Discussion

The major focus of our laboratory is the development and application of in vitro human-based experimental systems to allow accurate assessment of human drug properties. Our laboratory was one of the first to demonstrate successful cryopreservation of human hepatocytes (Loretz et al., 1989; Li, 2015b), and more recently human enterocytes (Ho et al., 2017), to retain viability and drug metabolism activities. We report here the characterization of a novel experimental system, namely, the MMHH experimental system—permeabilized, cofactor-supplemented cryopreserved human hepatocytes. MMHHs were developed in our laboratory to be an in vitro tool for drug metabolism studies with the desirable properties of both hepatocytes and subcellular fractions. As permeabilized, not lysed, cells, the MMHHs retain cellular integrity and subcellular structures and thereby contain complete enzyme pathways akin to CCHHs. Since the DME activities of MMHHs are driven by the cofactors—and therefore are not affected by viability—MMHHs, like subcellular fractions, are not subjected to handling damages, which can drastically affect the functions of CCHHs. The presence of all of the cellular organelles and the associated DMEs are the reasons that CCHHs are considered to be the gold standard in vitro drug metabolism system. These are also the attributes of MMHHs, but not subcellular fractions, which contain only a portion of the organelles, and therefore lack the DMEs associated with the missing organelles (Table 4).

A comparison of MMHHs and CCHHs in the activities of key DME pathways showed no apparent deficiencies in these activities in MMHHs. Both experimental systems were competent in phase 1 oxidation and phase 2 conjugation (including sequential events with oxidation followed by conjugation of the oxidative metabolites), as demonstrated using 7-ethoxycoumarin and acetaminophen as substrates. Similar cell density dependency

(from 0.25 to 2 million cells/ml), time course (up to 4 hours of incubation), and overall specific activities were observed for MMHHs and CCHHs in key phase 1 oxidation and phase 2 conjugation DME activities, with MMHHs competent in all pathways evaluated. Activities observed in MMHHs in general are within 2-fold of that for CCHHs (Fig. 7). Activities observed in MMHHs as percentage of that in CCHHs based on the rate of metabolite formation after an incubation duration of 30 minutes are: UDP-glucuronosyltransferase (271%), glutathione S-transferase (234%), SULT (204%), CYP1A2 (180%), CYP2A6 (177%), CYP2B6 (146%), CYP2C8 (131%), CYP2C9 (187%), CYP2E1 (132%), CYP3A4 (midazolam 1'-hydroxylation: (151%); testosterone 6 β -hydroxylation: (159%); CYP2J2 (100%); 7-ethoxycoumarin O-deethylation (179%); flavin-containing monooxygenase (169%); N-acetyltransferase 1 (134%); CYP2C19 (57%); CYP2D6 (86%); monoamine oxidase (89%); aldehyde oxidase (34%); carboxylesterase 2 (27%); acetaminophen glucuronidation (80%); and acetaminophen sulfation (43%).

Our results, therefore, show that MMHHs are active in all DME pathways evaluated with activities comparable to that for CCHHs, thereby suggesting that MMHHs can be used for drug metabolism studies that routinely employ CCHHs. The advantages of MMHHs over CCHHs include storage at -80°C instead of liquid nitrogen to retain viability and functions, direct use after thawing without a need for the time-consuming procedures of centrifugation cell recovery microscopy

TABLE 4
Organelle compositions of the various in vitro drug metabolism experimental systems

MMHHs, CCHHs, human liver microsomes (HLM) and postmitochondrial supernatant (PMS). The "X" indicates the presence of the organelles in the specified experimental systems. Organelle contents in various in vitro systems have been previously reviewed (Brandon et al., 2003).

| Organelle | Experimental System | | | |
|-----------------------|---------------------|------|-----|-----|
| | CCHH | MMHH | HLM | PMS |
| Plasma membranes | X | X | | |
| Cytosol | X | X | | X |
| Mitochondria | X | X | | |
| Lysosomes | X | X | | |
| Endoplasmic reticulum | X | X | X | X |
| Nucleus | X | X | | |

TABLE 5

Comparison of the procedures for drug metabolism studies with CCHHs and MMHHs

The time requirement for CCHHs will increase significantly if multiple individual lots (e.g., comparison of different hepatocyte donors) are used since each lot will require efforts in cell recovery, viability and cell concentration determination, and preparation of cell suspensions at the desired cell density. These additional steps are not required for MMHHs.

| Procedure | CCHH | MMHH |
|---|--|----------------------------------|
| Vial retrieval of from storage | From liquid nitrogen storage (3–5 minutes) | From –80°C storage (3–5 minutes) |
| Thawing of the cryopreserved hepatocytes in a 37°C water bath | 1.5–2 minutes | 1.5–2 minutes |
| Addition of the thawed hepatocytes to recovery medium followed by centrifugation for 10 minutes for cell recovery | 10–15 minutes | Not required |
| Microscopic examination of thawed hepatocytes for the quantification of viability and cell concentration | 10–15 minutes | Not required |
| Adjustment of cell suspension to contain 2× the final desired cell density | 5 minutes | Not required |
| Addition to vessels containing 2× test compounds ^a | | |
| Total time from freezer to incubation | 30.5–45 minutes | 7–10 minutes |

^aFor example, addition of 50 µl of hepatocytes to a 96-well plate with 50 µl of test articles at 2× the final concentration, followed by incubation for the desired time duration at 37°C: 1–3 minutes.

examination of cell viability, and cell concentration determination. We have estimated that the time required from thawing to incubation for drug metabolism evaluation is approximately 7–10 minutes for MMHHs versus 30.5–45 minutes for CCHHs (Table 5). Furthermore, unlike CCHHs in which cofactor compositions cannot be readily controlled, specific metabolic pathways can be evaluated with MMHHs based on the composition of the added cofactor mixtures. Selection of metabolic pathways based on cofactor compositions in MMHHs is demonstrated with coumarin metabolism where the use of UDPGA in the absence of PAPS led to formation of coumarin 7-glucuronide with no sulfate conjugate and vice versa.

Another potential advantage of MMHHs over CCHHs is that drug metabolism can be evaluated at cytotoxic drug concentrations. This was demonstrated with acetaminophen where drug metabolism activity is significantly decreased by drug cytotoxicity in CCHHs but not in MMHHs. This property of MMHHs should be useful in the evaluation of metabolic fate of highly toxic drugs such as anticancer therapeutics.

Our experimental findings suggest that MMHHs may represent a convenient hepatocyte-based experimental system that can be used for drug metabolism studies. Studies have been initiated in our laboratory on the application of MMHHs in the evaluation of metabolic stability and metabolite profiling, as well as reversible and mechanism-based DME inhibition. It needs to be emphasized, however, that CCHHs remain valuable and cannot be replaced by MMHHs in the evaluation of a number of key drug properties. These include uptake and efflux transport (Shitara et al., 2003; Bi et al., 2017; Izumi et al., 2017; Zhang et al., 2017), which require intact plasma membranes; enzyme induction (Fahmi et al., 2008; Zhou et al., 2011; Jones et al., 2017a), which requires RNA and protein synthesis; and in vitro hepatotoxicity (Zhang et al., 2016; Andersson, 2017; Proctor et al., 2017), which requires intact, viable hepatocytes. We believe that MMHHs should be complementary to current in vitro systems for the evaluation of human drug metabolism.

Acknowledgments

We gratefully acknowledge the technical assistance of Nick Ring in hepatocyte incubations and analytical chemistry sample preparation.

Authorship Contributions

Participated in research design: Li, Ho, Amaral, Loretz.

Conducted experiments: Ho, Amaral, Loretz.

Performed data analysis: Li, Ho.

Wrote or contributed to the writing of the manuscript: Li, Ho, Amaral, Loretz.

References

- Andersson TB (2017) Evolution of novel 3D culture systems for studies of human liver function and assessments of the hepatotoxicity of drugs and drug candidates. *Basic Clin Pharmacol Toxicol* **121**:234–238.
- Bi YA, Scialis RJ, Lazzaro S, Mathialagan S, Kimoto E, Keefer J, Zhang H, Vildhede AM, Costales C, Rodrigues AD, et al. (2017) Reliable rate measurements for active and passive hepatic uptake using plated human hepatocytes. *AAPS J* **19**:787–796.
- Brandon EF, Raap CD, Meijerman I, Beijnen JH, and Schellens JH (2003) An update on in vitro test methods in human hepatic drug biotransformation research: pros and cons. *Toxicol Appl Pharmacol* **189**:233–246.
- Chen S, Prieto Garcia L, Bergström F, Nordell P, and Grime K (2017) Intrinsic clearance assay incubational binding: a method comparison. *Drug Metab Dispos* **45**:342–345.
- Choi TA, Furimsky AM, Swezey R, Bunin DI, Byrge P, Iyer LV, Chang PY, and Aberger RJ (2015) In vitro metabolism and stability of the actinide chelating agent 3,4,3-L(1,2-HOPO). *J Pharm Sci* **104**:1832–1838.
- Di L, Atkinson K, Orozco CC, Funk C, Zhang H, McDonald TS, Tan B, Lin J, Chang C, and Obach RS (2013) In vitro-in vivo correlation for low-clearance compounds using hepatocyte relay method. *Drug Metab Dispos* **41**:2018–2023.
- Dinger J, Meyer MR, and Maurer HH (2014) Development of an in vitro cytochrome P450 cocktail inhibition assay for assessing the inhibition risk of drugs of abuse. *Toxicol Lett* **230**:28–35.
- Doshi U and Li AP (2011) Luciferin IPA-based higher throughput human hepatocyte screening assays for CYP3A4 inhibition and induction. *J Biomol Screen* **16**:903–909.
- Fabre G, Combalbert J, Berger Y, and Cano JP (1990) Human hepatocytes as a key in vitro model to improve preclinical drug development. *Eur J Drug Metab Pharmacokinet* **15**:165–171.
- Fahmi OA, Boldt S, Kish M, Obach RS, and Tremaine LM (2008) Prediction of drug-drug interactions from in vitro induction data: application of the relative induction score approach using cryopreserved human hepatocytes. *Drug Metab Dispos* **36**:1971–1974.
- Godoy P, Hewitt NJ, Albrecht U, Andersen ME, Ansari N, Bhattacharya S, Bode JG, Bolleyn J, Borner C, Böttger J, et al. (2013) Recent advances in 2D and 3D in vitro systems using primary hepatocytes, alternative hepatocyte sources and non-parenchymal liver cells and their use in investigating mechanisms of hepatotoxicity, cell signaling and ADME. *Arch Toxicol* **87**:1315–1530.
- Gombar VK, Silver IS, and Zhao Z (2003) Role of ADME characteristics in drug discovery and their in silico evaluation: in silico screening of chemicals for their metabolic stability. *Curr Top Med Chem* **3**:1205–1225.
- Halladay JS, Wong S, Jaffer SM, Sinhababu AK, and Khojasteh-Bakht SC (2007) Metabolic stability screen for drug discovery using cassette analysis and column switching. *Drug Metab Lett* **1**:67–72.
- Hewitt NJ, Lechón MJ, Houston JB, Hallifax D, Brown HS, Maurel P, Kenna JG, Gustavsson L, Lohmann C, Skonberg C, et al. (2007) Primary hepatocytes: current understanding of the regulation of metabolic enzymes and transporter proteins, and pharmaceutical practice for the use of hepatocytes in metabolism, enzyme induction, transporter, clearance, and hepatotoxicity studies. *Drug Metab Rev* **39**:159–234.
- Ho MD, Ring N, Amaral K, Doshi U, and Li AP (2017) Human enterocytes as an in vitro model for the evaluation of intestinal drug metabolism: characterization of drug-metabolizing enzyme activities of cryopreserved human enterocytes from twenty-four donors. *Drug Metab Dispos* **45**:686–691.
- Iwatsubo T, Suzuki H, Shimada N, Chiba K, Ishizaki T, Green CE, Tyson CA, Yokoi T, Kamataki T, and Sugiyama Y (1997) Prediction of in vivo hepatic metabolic clearance of YM796 from in vitro data by use of human liver microsomes and recombinant P-450 isozymes. *J Pharmacol Exp Ther* **282**:909–919.
- Izumi S, Nozaki Y, Komori T, Takenaka O, Maeda K, Kusuhashi H, and Sugiyama Y (2017) Comparison of the predictability of human hepatic clearance for organic anion transporting polypeptide substrate drugs between different in vitro-in vivo extrapolation approaches. *J Pharm Sci* **106**:2678–2687.
- Jones BC, Rollison H, Johansson S, Kanebratt KP, Lambert C, Vishwanathan K, and Andersson TB (2017a) Managing the risk of CYP3A induction in drug development: a strategic approach. *Drug Metab Dispos* **45**:35–41.
- Jones BC, Srivastava A, Colclough N, Wilson J, Reddy VP, Amberntsson S, and Li D (2017b) An investigation into the prediction of in vivo clearance for a range of flavin-containing monooxygenase substrates. *Drug Metab Dispos* **45**:1060–1067.
- Joo J, Lee B, Lee T, and Liu KH (2014) Screening of six UGT enzyme activities in human liver microsomes using liquid chromatography/triple quadrupole mass spectrometry. *Rapid Commun Mass Spectrom* **28**:2405–2414.

- Kazmi F, Barbara JE, Yerino P, and Parkinson A (2015) A long-standing mystery solved: the formation of 3-hydroxydesloratadine is catalyzed by CYP2C8 but prior glucuronidation of desloratadine by UDP-glucuronosyltransferase 2B10 is an obligatory requirement. *Drug Metab Dispos* **43**:523–533.
- Li AP (2007) Human hepatocytes: isolation, cryopreservation and applications in drug development. *Chem Biol Interact* **168**:16–29.
- Li AP (2015a) inventor, Li AP assignee. Cell preparation method. US patent no. U.S. patent 9,078,430 B2. 2015 Jul 14.
- Li AP (2015b) Evaluation of adverse drug properties with cryopreserved human hepatocytes and the integrated discrete multiple organ co-culture (IdMOC™) system. *Toxicol Res* **31**: 137–149.
- Li AP and Doshi U (2011) Higher throughput human hepatocyte assays for the evaluation of time-dependent inhibition of CYP3A4. *Drug Metab Lett* **5**:183–191.
- Loretz LJ, Li AP, Flye MW, and Wilson AG (1989). Optimization of cryopreservation procedures for rat and human hepatocytes. *Xenobiotica* **19**:489–498.
- Obach RS (2011) Predicting clearance in humans from in vitro data. *Curr Top Med Chem* **11**: 334–339.
- Peng CC, Doshi U, Prakash C, and Li AP (2016) A novel plated hepatocyte relay assay (PHRA) for in vitro evaluation of hepatic metabolic clearance of slowly metabolized compounds. *Drug Metab Lett* **10**:3–15.
- Proctor WR, Foster AJ, Vogt J, Summers C, Middleton B, Pilling MA, Shienson D, Kijanska M, Ströbel S, Kelm JM, et al. (2017) Utility of spherical human liver microtissues for prediction of clinical drug-induced liver injury. *Arch Toxicol* **91**:2849–2863.
- Shitara Y, Li AP, Kato Y, Lu C, Ito K, Itoh T, and Sugiyama Y (2003) Function of uptake transporters for taurocholate and estradiol 17 β -D-glucuronide in cryopreserved human hepatocytes. *Drug Metab Pharmacokinet* **18**:33–41.
- Ulrich RG, Bacon JA, Cramer CT, Peng GW, Petrella DK, Stryd RP, and Sun EL (1995) Cultured hepatocytes as investigational models for hepatic toxicity: practical applications in drug discovery and development. *Toxicol Lett* **82–83**:107–115.
- Walsky RL, Bauman JN, Bourcier K, Giddens G, Lapham K, Negahban A, Ryder TF, Obach RS, Hyland R, and Goosen TC (2012) Optimized assays for human UDP-glucuronosyltransferase (UGT) activities: altered alamethicin concentration and utility to screen for UGT inhibitors. *Drug Metab Dispos* **40**:1051–1065.
- Wood FL, Houston JB, and Hallifax D (2017) Clearance prediction methodology needs fundamental improvement: trends common to rat and human hepatocytes/microsomes and implications for experimental methodology. *Drug Metab Dispos* **45**:1178–1188.
- Zhang J, Doshi U, Suzuki A, Chang CW, Borlak J, Li AP, and Tong W (2016) Evaluation of multiple mechanism-based toxicity endpoints in primary cultured human hepatocytes for the identification of drugs with clinical hepatotoxicity: results from 152 marketed drugs with known liver injury profiles. *Chem Biol Interact* **255**:3–11.
- Zhang Y, Jackson JP, St Claire RL, III, Freeman K, Brouwer KR, and Edwards JE (2017) Obeticholic acid, a selective farnesoid X receptor agonist, regulates bile acid homeostasis in sandwich-cultured human hepatocytes. *Pharmacol Res Perspect* **5**:e00329.
- Zhou D, Andersson TB, and Grimm SW (2011) In vitro evaluation of potential drug-drug interactions with Ticagrelor: cytochrome P450 reaction phenotyping, inhibition, induction, and differential kinetics. *Drug Metab Dispos* **39**:703–710.

Address correspondence to: Dr. Albert P. Li, 9221 Rumsey Road Suite 8, Columbia, MD 21045. E-mail: lialbert@invitroadmet.com
