

FUJIFILM**Wako**

Code No. 207-19982 (25 g)
209-19981 (100 g)
201-19985 (500 g)

Trypsin 1 : 250, from Porcine Pancreas

トリプシン 1 : 250, プタ膵臓由来

Trypsin is a serine protease that cleaves peptides at the C-terminus of basic amino acid residues (Lysine and Arginine).

This product is for laboratory research use only ; use in any such application is the responsibility of the user.

Source : Porcine Pancreas


Trypsin Activity : Indicated on the label

CAS No. : 9002-07-7

EC No. : 3.4.21.4

Assay Methods

1. Principle

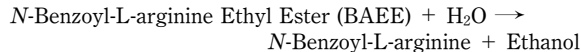

2. Reagents

A. 0.067 mol/L Phosphate Buffer, pH 7.6

B. BAEE Solution

Dissolve 86 mg of *N*-Benzoyl-L-arginine Ethyl Ester Hydrochloride in 100 mL of water. Transfer 10 mL of the solution and dilute with Reagent A to bring the volume up to 100 mL.

C. 0.001 mol/L HCl

D. Trypsin Solution

Dissolve 50 mg of this product in 100 mL of Reagent C. Transfer 2.5 mL of the solution and dilute with Reagent C to bring the volume up to 100 mL.

3. Procedures

Reagent	Test	Blank
B	3.0 mL	3.0 mL
Pre-incubate at 25 °C for 5 minutes		
D	0.2 mL	—
C	—	0.2 mL
Immediately, measure the absorbance of 253 nm at 25 °C for 5 minutes with water as reference.		

— 1/2 —

4. Unit Definition

1 USP Unit is the amount of enzyme which increases 0.003 of absorbance per minute under the above conditions.

(Calculation)

$$A = (E_1 - E_2) \times \frac{1}{0.003} \times \frac{100}{0.2 \times S} \times \frac{100}{2.5}$$

A : Trypsin Activity (USP Units/mg)

E_1 : Absorbance change per minute in the test

E_2 : Absorbance change per minute in the blank

0.003 : Absorbance change per minute per 1 USP Unit

0.2 : The amount of diluted trypsin solution (mL)

S : The amount of weighed this product (mg)

Trans *

4. 単位の定義

上記の条件で、1 分間に吸光度を 0.003 増加させる酵素量を 1 USP Unit とする。

(計 算)

$$A = (E_1 - E_2) \times \frac{1}{0.003} \times \frac{100}{0.2 \times S} \times \frac{100}{2.5}$$

A : トリプシン活性 (USP Units/mg)

E_1 : 本試験の 1 分間あたりの吸光度変化

E_2 : 空試験の 1 分間あたりの吸光度変化

0.003 : 1 USP Unit あたりの吸光度変化/分

0.2 : 本品溶液の使用量 (mL)

S : 量りとした本品量 (mg)

Storage : Store at room temperature

Package :

Code No.	Packaging
207-19982	25 g
209-19981	100 g
201-19985	500 g

* : **Trans** is the Japanese translation.

FUJIFILM Wako Pure Chemical Corporation

1-2, Doshomachi 3-Chome, Chuo-Ku, Osaka 540-8605, Japan

Telephone : + 81-6-6203-3741

Facsimile : + 81-6-6201-5964

<http://fwk.fujifilm.co.jp>

FUJIFILM Wako Chemicals U.S.A. Corporation

1600 Bellwood Road

Richmond, VA 23237

U.S.A.

Telephone : + 1-804-271-7677

Facsimile : + 1-804-271-7791

<http://www.wakousa.com>

FUJIFILM Wako Chemicals Europe GmbH

Fuggenstrasse 12

D-41468 Neuss

Germany

Telephone : + 49-2131-3111-0

Facsimile : + 49-2131-311100

<http://www.wako-chemicals.de>

2304KA2

— 2/2 —